

Instructions for use

Pneumocystis ELITe MGB® Kit

reagents for DNA Real-Time PCR



REF RTS150ING

UDI 08033891486723



CHANGE HISTORY

| Revision | Notice of change | Date (dd/mm/yyyy) |
|----------|---|-------------------|
| 02 | Extended use of the product in association with «ELiTe BeGenius®» instrument (REF INT040) and Sputum matrix. Update of PERFORMANCE CHARACTERISTICS: — update of ULoQ value for sputum matrix — update of Linear measuring range for sputum performed in matrix instead of PBS. New graphics and content setting of the IFU. | 25/07/2024 |
| 01 | Extended use of the product on «ELiTe BeGenius®» instrument (REF INT040) and BAL matrix. Updated ULoQ/LLoQ value calculated on BAL matrix. Description of IC cut off value. | 15/05/2023 |
| 00 | New product development | 09/11/2019 |

NOTE

The revision of this IFU is also compatible with the previous versions of the kit

TABLE OF CONTENT

| | |
|--|-----------|
| 1 INTENDED USE | 4 |
| 2 ASSAY PRINCIPLE | 4 |
| 3 PRODUCT DESCRIPTION | 4 |
| 4 MATERIALS PROVIDED IN THE PRODUCT | 4 |
| 5 MATERIALS REQUIRED BUT NOT PROVIDED IN THE PRODUCT..... | 5 |
| 6 OTHER PRODUCTS REQUIRED..... | 5 |
| 7 WARNINGS AND PRECAUTIONS | 5 |
| 8 SPECIMENS AND CONTROLS | 7 |
| 9 ELiTe InGenius PROCEDURE..... | 9 |
| 10 ELiTe BeGenius PROCEDURE | 15 |
| 11 PERFORMANCE CHARACTERISTICS..... | 19 |
| 12 REFERENCES..... | 28 |
| 13 PROCEDURE LIMITATIONS | 28 |
| 14 TROUBLESHOOTING | 29 |
| 15 SYMBOLS | 31 |
| 16 NOTICE TO PURCHASER: LIMITED LICENSE | 32 |
| Appendix A QUICK START GUIDE..... | 33 |

1 INTENDED USE

The product **Pneumocystis ELiTe MGB Kit** is an *in vitro* diagnostic medical device intended to be used by healthcare professionals as a quantitative nucleic acids Real-Time PCR assay for the **detection and the quantification of the DNA of *Pneumocystis jirovecii* (PJ)** in DNA samples extracted from clinical specimens.

The assay is validated in association with the **ELiTe InGenius®** and **ELiTe BeGenius®** instruments, automated and integrated systems for extraction, Real-Time PCR and results interpretation, using human Broncho-alveolar lavage (BAL) / bronchial aspirate (BA) and liquefied Sputum specimens

The product is intended for *in vitro* diagnostic use as an aid in diagnosis and monitoring of *Pneumocystis jirovecii* infection.

The results must be interpreted in combination with all relevant clinical observations and laboratory outcomes.

2 ASSAY PRINCIPLE

The assay is a quantitative Real-Time PCR detecting *Pneumocystis jirovecii* (**PJ**) DNA, isolated from specimens and amplified using the assay reagent **PJ PCR Mix** that contains primers and probes with ELiTe MGB® technology.

The ELiTe MGB probes are activated when hybridize with the related PCR products. **ELiTe InGenius** and **ELiTe BeGenius** monitor fluorescence increase and calculate the threshold cycles (Ct) and the melting temperatures (Tm). The PJ DNA quantity is calculated based on a stored calibration curve.

In the ELiTe MGB probes the fluorophores are quenched in the random-coiled, single-stranded state of probe. The fluorophores are active in the probe / amplicon duplex as the quencher is spatially separated from the fluorophore. Note the fluorophore is not cleaved during PCR and can be utilized for dissociation analysis and melting temperature calculation.

3 PRODUCT DESCRIPTION

The **Pneumocystis ELiTe MGB Kit** provides the assay reagent **PJ PCR Mix**, an optimized and stabilized PCR mixture that contains the specific primers and probes for:

- The **mitochondrial Large Subunit of rRNA gene (mtLSU)** of PJ, detected in Channel **PJ**; the probe is stabilized by MGB, quenched by the Eclipse Dark Quencher®, and labelled by FAM dye,
- Internal Control, artificial sequence **IC2**, detected in Channel **IC**; the probe is stabilized by MGB, quenched by the Eclipse Dark Quencher, and labelled by AquaPhluor 525 (AP525) dye.

The **PJ PCR Mix** also contains buffer, magnesium chloride, triphosphate nucleotides and hot start DNA Polymerase.

The **Pneumocystis ELiTe MGB Kit** contains sufficient reagents for 96 tests in association with **ELiTe InGenius** and **ELiTe BeGenius (12 tests each tube)**, with 20 µL used per reaction.

The **Pneumocystis ELiTe MGB Kit** can be also used in association with other equivalent instruments.

4 MATERIALS PROVIDED IN THE PRODUCT

Table 1

| Component | Description | Quantity | Classification of hazards |
|-------------------------------------|--|-------------------|---------------------------|
| PJ PCR Mix ref. RTS150ING | Mixture of reagents for Real-Time PCR in tube with WHITE cap | 8 x 280 µL | - |

5 MATERIALS REQUIRED BUT NOT PROVIDED IN THE PRODUCT

- Laminar airflow hood.
- Disposable nitrile powder-free gloves or similar material.
- Vortex mixer.
- Bench centrifuge (~5,000 RPM).
- Bench microcentrifuge (~13,000 RPM).
- Micropipettes and sterile tips with aerosol filter or sterile positive displacement tips (0.5-10 µL, 2-20 µL, 5-50 µL, 50-200 µL, 200-1000 µL).
- 2.0 mL sterile screw capped tubes (Sarstedt, Germany, ref. 72.694.005).
- Molecular biology grade water.

6 OTHER PRODUCTS REQUIRED

The reagents for the extraction of sample DNA, the extraction and inhibition internal control, the amplification positive and negative controls, the DNA standards and the consumables **are not** provided with this product.

For automated extraction of nucleic acids, Real-Time PCR and result interpretation of samples, the following products are required:

Table 2

| Instruments and Softwares | Products and Reagents |
|--|--|
| ELiTe InGenius (ELiTechGroup S.p.A., EG SpA ref. INT030) ELiTe InGenius Software version 1.3.0.19 (or later) PJ ELiTe _STD , Assay Protocol with parameters for Calibrators analysis PJ ELiTe _PC , Assay Protocol with parameters for Positive Control analysis PJ ELiTe _NC , Assay Protocol with parameters for Negative Control analysis PJ ELiTe _BAL_200_100 , Assay Protocol with parameters for BAL/BA specimen analysis PJ ELiTe _SP_200_100 , Assay Protocol with parameters for sputum specimen analysis | ELiTe InGenius SP200 (EG SpA, ref. INT032SP200) ELiTe InGenius SP 200 Consumable Set (EG SpA, ref. INT032CS) ELiTe InGenius PCR Cassette (EG SpA, ref. INT035PCR), ELiTe InGenius Waste Box (EG SpA, ref. F2102-000) 300 µL Filter Tips Axygen (Corning Life Sciences Inc., ref. TF-350-L-R-S) with ELiTe InGenius only 1000 µL Filter Tips Tecan (Tecan, Switzerland, ref. 30180118) with ELiTe BeGenius only Pneumocystis - ELiTe Positive Control (EG SpA, ref. CTR150PLD). Pneumocystis ELiTe Standard (EG SpA, ref. STD150PLD) CPE - Internal Control (EG SpA, ref. CTRCPE). |
| ELiTe BeGenius (EG SpA ref. INT040) ELiTe BeGenius Software version 2.2.1. (or later) PJ ELiTe _Be_STD , Assay Protocol with parameters for Calibrators analysis PJ ELiTe _Be_PC , Assay Protocol with parameters for Positive Control analysis PJ ELiTe _Be_NC , Assay Protocol with parameters for Negative Control analysis PJ ELiTe _Be_BAL_200_100 , Assay Protocol with parameters for BAL/BA specimen analysis PJ ELiTe _Be_SP_200_100 , Assay Protocol with parameters for sputum specimen analysis | |

7 WARNINGS AND PRECAUTIONS

This product is designed for in-vitro use only.

7.1 General warnings and precautions

Handle and dispose of all biological samples as if they were infectious. Avoid direct contact with biological samples. Avoid splashing or spraying. Tubes, tips and other materials that come into contact with the biological samples must be treated for at least 30 minutes with 3% sodium hypochlorite (bleach) or autoclaved for one hour at 121°C before disposal.

Handle and dispose of all reagents and all materials used to carry out the assay as if they were infectious. Avoid direct contact with the reagents. Avoid splashing or spraying. Waste must be handled and disposed of in compliance with adequate safety standards. Disposable combustible material must be incinerated. Liquid waste containing acids or bases must be neutralized before disposal. Do not allow extraction reagents to contact sodium hypochlorite (bleach).

- Wear suitable protective clothes and gloves and protect eyes and face.
- Never pipette solutions by mouth.
- Do not eat, drink, smoke or apply cosmetic products in the work areas.
- Carefully wash hands after handling samples and reagents.
- Dispose of leftover reagents and waste in compliance with the regulations in force.
- Carefully read all the instructions provided before running the assay.
- While running the assay, follow the product instructions provided.
- Do not use the product after the indicated expiry date.
- Only use reagents provided with the product and those recommended by the manufacturer.
- Do not use reagents from different batches.
- Do not use reagents from other manufacturers.

7.2 Warnings and precautions for molecular biology

Molecular biology procedures require qualified and trained staff to avoid the risk of erroneous results, especially due to sample nucleic acid degradation or sample contamination by PCR products.

Never transfer lab coats, gloves or tools from the area designated for the amplification / detection of amplification products to the area designated for the extraction / preparation of the amplification reactions.

Laboratory coats, gloves and tools dedicated to work session setup are needed.

The samples must be suitable and, if possible, dedicated for this type of analysis. Samples must be handled under a laminar airflow hood. Pipettes used to handle samples must be exclusively used for this specific purpose. The pipettes must be of the positive displacement type or be used with aerosol filter tips. The tips used must be sterile, free from DNases and RNases, and free from DNA and RNA.

The reagents must be handled under a laminar airflow hood. The pipettes used to handle the reagents must be exclusively used for this purpose. The pipettes must be of the positive displacement type or be used with aerosol filter tips. The tips used must be sterile, free from DNases and RNases, and free from DNA and RNA.

The extraction products must be handled in such a way as to minimize dispersion into the environment in order to avoid the possibility of contamination.

The PCR Cassette must be handled carefully and never opened to avoid PCR product diffusion into the environment and sample and reagent contamination.

7.3 Warnings and precautions specific for the components

Table 3

| Component | Storage temperature | Use from first opening | Freeze / Thaw cycles | On board stability (ELiTe InGenius and ELiTe BeGenius) |
|-------------|---|------------------------|----------------------|---|
| PJ— PCR Mix | -20 °C or below (protected from light) | one month | up to seven | up to seven separate* sessions of three hours each or up to 7 consecutive hours (2 sessions of 3 hours each and the time needed to start a third session) |

* with intermediate freezing

8 SPECIMENS AND CONTROLS

8.1 Specimens

This product is intended for use on the **ELiTe InGenius** and **ELiTe BeGenius** with the following clinical specimens identified and handled according to laboratory guidelines, and collected, transported, and stored under the following conditions:

Table 4

| Specimen | Collection requirements | Transport/Storage conditions | | |
|--|---|------------------------------|-------------|-------------|
| | | +2 / +8 °C | -20 ± 10 °C | -70 ± 15 °C |
| Bronchoalveolar Lavage / bronchial aspirate (BAL/BA) | in sterile physiological solution or sterile PBS* | ≤1 week | ≤ 30 days | ≤ 1 year |
| Sputum (SP) | | ≤1 week | ≤ 30 days | long period |

*PBS: Phosphate Buffered Saline solution

If BAL/BA samples are particularly mucous, they can be liquefied by dithiothreitol based reagents (e.g. Sputasol, Oxoid, Thermo Fisher Scientific) as per laboratory guidelines.

Samples of Sputum must be liquefied by dithiothreitol based reagents (e.g. Sputasol, Oxoid, Thermo Fisher Scientific) as per laboratory guidelines.

It is recommended to divide the specimens into aliquots before freezing to prevent repeated freeze / thaw cycles. When using frozen samples, thaw the samples just before the extraction to avoid possible nucleic acid degradation.

To perform samples testing on the **ELiTe InGenius** and **ELiTe BeGenius**, the following Assay Protocols must be used. These IVD protocols were specifically validated with ELiTe MGB Kits and the **ELiTe InGenius** or **ELiTe BeGenius** with the indicated matrices.

Table 5 Assay Protocols for Pneumocystis ELiTe MGB Kit

| Specimen | Instrument | Assay Protocol Name | Report | Characteristics |
|--|----------------|-------------------------|-----------|--|
| Bronchoalveolar Lavage / bronchial aspirate (BAL/BA) | ELiTe InGenius | PJ ELiTe_BAL_200_100 | copies/mL | Extraction Input Volume: 200 µL Extraction Elution Volume: 100 µL Internal Control: 10 µL Sonication: NO Dilution Factor: 1 PCR Mix volume: 20 µL Sample PCR input volume: 20 µL |
| | ELiTe BeGenius | PJ ELiTe_Be_BAL_200_100 | | |
| Sputum (SP) | ELiTe InGenius | PJ ELiTe_SP_200_100 | copies/mL | Extraction Input Volume: 200 µL Extraction Elution Volume: 100 µL Internal Control: 10 µL Sonication: NO Dilution Factor: 1 PCR Mix volume: 20 µL Sample PCR input volume: 20 µL |
| | ELiTe BeGenius | PJ ELiTe_Be_SP_200_100 | | |

For all protocols, 200 µL of sample must be transferred into Extraction tube (for ELiTe InGenius) or 2 mL Sarstedt Tube (for ELiTe BeGenius).

NOTE

Pipetting samples to the **Extraction tube** or to the **2 mL Sarstedt Tube** might **generate contamination**. Use the appropriate pipettes and follow all recommendations reported in the “Warnings and Precautions” section.

Purified nucleic acids can be left at room temperature for 16 hours and stored at -20 °C or below for no longer than one month.

Refer to “Potentially Interfering Substances” in the Performance Characteristics section to check data concerning interfering substances.

Quantities of human genomic DNA higher than 1 µg extracted from the sample could inhibit the real-time amplification

8.2 PCR calibrators and controls

The Calibration curve must be generated and approved for each lot of PCR reagent.

- For the calibration curve, use the four levels of the product **Pneumocystis ELiTe Standard** (not provided with this kit) with the **PJ ELiTe_STD** or **PJ ELiTe_Be_STD** Assay Protocols.

NOTE

The concentration of Q – PCR Standards are expressed in copies / reaction (10⁵ copies / rxn, 10⁴ copies / rxn, 10³ copies / rxn, 10² copies / rxn).

PCR control results must be generated and approved for each lot of PCR reagent.

- For the Positive Control, use the product **Pneumocystis - ELiTe Positive Control** (not provided with this kit) with the **PJ ELiTe_PC** or **PJ ELiTe_Be_PC** Assay Protocols,
- For the Negative Control, use molecular biology grade water (not provided with this kit) with the **PJ ELiTe_NC** or **PJ ELiTe_Be_NC** Assay Protocols.

NOTE

The **ELiTe InGenius** and **ELiTe BeGenius** allow generation and storage of the calibration curve and PCR control validation for each lot of PCR reagent.

Calibration curves expire after **60 days**, at which time it is necessary to re-run the calibration.

PCR control results expire after **15 days**, at which time it is necessary to re-run the Positive and Negative Controls.

The Calibrators and PCR controls must be re-run if any of the following events occur:

- a new lot of reagents is used,
- results of quality control analysis (see following paragraph) are out of specification,
- any major maintenance or service is performed on the **ELiTe InGenius** or **ELiTe BeGenius** instruments.

8.3 Quality controls

Verification of the extraction and PCR procedure is recommended. Archived samples or certified reference material may be used. External controls should be used in accordance with local, state, and federal accrediting organizations, as applicable.

9 ELiTe InGenius PROCEDURE

The procedure to use the **Pneumocystis ELiTe MGB Kit** with the **ELiTe InGenius** consists of three steps:

Table 6

| | | |
|--------|--------------------------------------|--|
| STEP 1 | Verification of the system readiness | |
| STEP 2 | Session setup | A) Sample run (Extract + PCR) |
| | | B) Eluted sample run (PCR Only) |
| | | C) Calibration run (PCR Only) |
| | | D) Positive Control and Negative Control run (PCR Only) |
| STEP 3 | Review and approval of results | A) Validation of Calibration curve |
| | | B) Validation of Positive Control and Negative Control results |
| | | C) Validation of sample results |
| | | D) Sample result reporting |

9.1 STEP 1 – Verification of the system readiness

Before starting the session:

- switch on the **ELiTe InGenius** and login in “**CLOSED**” mode,
- in the “Calibration” menu on the Home page, verify the Calibrators (**Q - PCR Standard**) are approved and valid (Status) for the **PCR Mix** lot to be used. If no valid Calibrators are available for the **PCR Mix** lot, perform calibration as described in the following sections,
- in the “Controls” menu on the Home page, verify that the PCR Controls (**Positive Control, Negative Control**) are approved and valid (Status) for the **PCR Mix** lot to be used. If no valid PCR Controls are available for the **PCR Mix** lot, run the PCR Controls as described in the following sections,
- choose the type of run, following the instructions on the Graphical User Interface (GUI) for the session setup and using the Assay Protocols provided by EG SpA (see “[Specimens and Controls](#)”).

If the Assay Protocol of interest is not loaded in the system, contact your local ELiTechGroup Customer Service.

9.2 STEP 2 – Session Setup

The **Pneumocystis ELiTe MGB Kit** can be used on **ELiTe InGenius** to perform:

- A. Sample run (Extract + PCR),
- B. Eluted sample run, (PCR Only),
- C. Calibration run (PCR Only),
- D. Positive Control and Negative Control run (PCR Only).

All required parameters are included in the Assay Protocols available on the instrument and are loaded automatically when the Assay Protocol is selected.

NOTE

The **ELiTe InGenius** can be connected to the “Laboratory Information System” (LIS) which enables downloading the session information. Refer to the instrument manual for more details.

Before to setup a run:

Thaw the needed **PCR Mix** tubes at room temperature for 30 minutes. Each tube is sufficient for **12 tests** in optimized conditions (2 or more tests per session). Mix gently, then spin down the contents for 5 seconds and keep on ice or cool block.

NOTE

Protect the **PCR Mix** from light while thawing because this reagent is photosensitive.

To set up one of the four types of run follow the steps below while referring to the GUI:

| | A. Sample run (Extract + PCR) | B. Eluted sample run (PCR Only) |
|---|--|--|
| 1 | Identify samples and, if needed, thaw at room temperature, mix gently, spin down the contents for 5 seconds and keep on ice or cool block. If required, transfer 200 µL of sample in an Extraction tube previously labelled. Thaw the needed CPE tubes at room temperature for 30 minutes. Mix gently, spin down the contents for 5 seconds and keep on ice or cool block. Each tube is sufficient for 12 extractions. | Thaw the Elution tube containing the extracted nucleic acids at room temperature. Mix gently, then spin down the contents for 5 seconds and keep on ice or cool block. |
| 2 | Select “ Perform Run ” from the “Home” screen. | Select “ Perform Run ” from the “Home” screen. |
| 3 | Ensure the “Extraction Input Volume” is 200 µL and the “Extracted Elute Volume” is 100 µL. | Ensure the “Extraction Input Volume” is 200 µL and the “Extracted Elute Volume” is 100 µL. |
| 4 | For each sample, assign a Track and enter the “SampleID” (SID) by typing or by scanning the sample barcode. | For each sample, assign a Track and enter the “SampleID” (SID) by typing or by scanning the sample barcode. |
| 5 | Select the Assay Protocol in the “Assay” column (see “Specimens and Controls”). | Select the Assay Protocol in the “Assay” column (see “Specimens and Controls”). |
| 6 | Ensure the “Protocol” displayed is: “Extract + PCR”. | Select “PCR Only” in the “Protocol” column. |
| 7 | Select the sample loading position as “Extraction Tube” in the “Sample Position” column. Ensure the “ Dilution factor ” is “1”. | Ensure the sample loading position in the “Sample Position” column is “Elution Tube (bottom row)”. Ensure the “ Dilution factor ” is “1”. |
| 8 | Click “Next” to continue. | Click “Next” to continue. |
| 9 | Load CPE and the PCR Mix on the “Inventory Block” referring to the “Load List” and enter CPE and PCR Mix lot number, expiry date and number of reactions for each tube. | Load the PCR Mix on the “Inventory Block” referring to the “Load List” and enter PCR Mix lot number, expiry date and number of reactions for each tube. |

| | A. Sample run (Extract + PCR) | B. Eluted sample run (PCR Only) |
|----|--|--|
| 10 | Click "Next" to continue. | Click "Next" to continue. |
| 11 | Verify the tips in the "Tip Racks" in the "Inventory Area" and replace Tip Racks if necessary. | Verify the tips in the "Tip Racks" in the "Inventory Area" and replace Tip Racks if necessary. |
| 12 | Click "Next" to continue. | Click "Next" to continue. |
| 13 | Load PCR Cassette, ELITE InGenius SP 200 extraction cartridges, and all required consumables and samples to be extracted. | Load PCR Cassette and Elution tubes with samples extracted. |
| 14 | Click "Next" to continue. | Click "Next" to continue. |
| 15 | Close the instrument door. | Close the instrument door. |
| 16 | Press "Start". | Press "Start". |

| | C. Calibration run (PCR Only) | D. Positive Control and Negative Control run (PCR Only) |
|----|--|--|
| 1 | Thaw the needed Q-PCR Standard tubes (Cal1: Q-PCR Standard 10 ² , Cal2: Q-PCR Standard 10 ³ , Cal3: Q-PCR Standard 10 ⁴ , Cal4: Q-PCR Standard 10 ⁵) at room temperature for 30 minutes. Mix gently, then spin down the contents for 5 seconds and keep on ice or cool block. | Thaw Positive Control tubes at room temperature for 30 minutes. Mix gently, then spin down the contents for 5 seconds and keep on ice or cool block. Prepare the Negative Control by transferring at least 50 µL of molecular biology grade water to an "Elution tube", provided with the ELITE InGenius SP 200 Consumable Set. |
| 2 | Select "Perform Run" from the "Home" screen. | Select "Perform Run" from the "Home" screen. |
| 3 | Ensure the "Extraction Input Volume" is 200 µL and the "Extracted Elute Volume" is 100 µL. | Ensure "Extraction Input Volume" is 200 µL and "Extracted Elute Volume" is 100 µL. |
| 4 | For the Q-PCR Standard, assign the "Track", select the Assay Protocol (see "Specimen and Controls") in the "Assay" column and enter the reagent lot number and expiry date. | Select the Assay Protocol in the "Assay" column (see "Specimens and Controls"). Enter the lot number and expiry date of the Positive Control and of the molecular biology grade water. |
| 5 | Ensure "PCR Only" is selected in the "Protocol" column. | Ensure "PCR Only" is selected in the "Protocol" column. |
| 6 | Ensure the sample loading position in "Sample Position" column is "Elution Tube (bottom row)". | Ensure the sample loading position in the "Sample Position" column is "Elution Tube (bottom row)". |
| 7 | Load the PCR Mix on the "Inventory Block" referring to the Load List and enter the PCR Mix lot number, expiry date and number of reactions for each tube. | Load the PCR Mix on the "Inventory Block" referring to the "Load List" and enter the PCR Mix lot number, expiry date and number of reactions for each tube. |
| 8 | Click "Next" to continue. | Click "Next" to continue. |
| 9 | Verify the tips in the "Tip Rack(s)" in the "Inventory Area" and replace Tip Rack(s) if necessary. | Verify the tips in the "Tip Rack(s)" in the "Inventory Area" and replace Tip Rack(s) if necessary. |
| 10 | Click "Next" to continue. | Click "Next" to continue. |
| 11 | Load the PCR Cassette and the Q-PCR Standard tubes. | Load PCR Cassette, Positive Control and Negative Control. |
| 12 | Click "Next" to continue. | Click "Next" to continue. |
| 13 | Close the instrument door. | Close the instrument door. |
| 14 | Press "Start" | Press "Start". |

When the session is finished, the **ELITE InGenius** allows users to view, approve, store the results, print and save the report.

NOTE

At the end of the run the remaining Extracted Sample in the **Elution tube** must be removed from the instrument, capped, identified and stored at -20 ± 10 °C for no longer than one month. Avoid spilling of the Extracted Sample.

NOTE

At the end of the run the **PCR Mix** can be removed from the instrument, capped and stored at -20 °C or below or can be kept on board in the refrigerated block for up to 7 hours (2 sessions of 3 hours each and the time needed to start a third session); mix gently and spin down the content for 5 seconds before starting the next session.

NOTE

At the end of the run the remaining **Q - PCR Standard** can be removed from the instrument, capped, and stored at -20 °C or below. Avoid spilling the Q - PCR Standard.

NOTE

The **Q - PCR Standard** can be used for 4 separate sessions of 2 hours each.

NOTE

At the end of the run the remaining **Positive Control** can be removed from the instrument, capped and stored at -20 °C or below. Avoid the spilling of the Positive Control. The remaining **Negative Control** must be discarded.

NOTE

The **Positive Control** can be used for 4 separate sessions of 3 hours each.

NOTE

At the end of the run the **PCR Cassette** and the other consumables must be disposed of following all governmental and environmental regulations. Avoid spilling the reaction products.

9.3 STEP 3 - Review and approval of results

The **ELITE InGenius** monitors target and internal control fluorescence signals for each reaction and automatically applies the Assay Protocol parameters to generate PCR curves which are then interpreted into results.

At the end of the run, the “Results Display” screen is automatically shown. In this screen the results and the run information are shown. From this screen results can be approved, and reports printed or saved (“Sample Report” or “Track Report”). Refer to the instrument manual for more details.

NOTE

The **ELITE InGenius** can be connected to the “Laboratory Information System” (LIS) which enables uploading the session results to the laboratory data center. Refer to the instrument manual for more details.

The **ELITE InGenius** generates results with the **Pneumocystis ELITE MGB Kit** through the following procedure:

- A. Validation of Calibration curve,
- B. Validation of Positive Control and Negative Control results,
- C. Validation of sample results,
- D. Sample result reporting.

9.3.1 A. Validation of Calibration curve

The **ELiTe InGenius software** interprets the PCR results for the target of the Calibrator reactions with the **ELiTe STD** Assay Protocol parameters. The resulting Ct versus concentration produces the Calibration curve.

The Calibration curves, specific for the PCR reagent lot, are recorded in the database (Calibration). They can be viewed and approved by "Administrator" or "Analyst" users, following the GUI instructions.

The Calibration curve expires **after 60 days**.

NOTE

If the Calibration curve does not meet the acceptance criteria, the "Failed" message is shown on the "Calibration" screen. In this case, the results cannot be approved and the Calibrator amplification reactions must be repeated. In addition, if samples were included in the run, these are not quantified and must also be repeated to generate quantitative results.

9.3.2 B. Validation of amplification Positive Control and Negative Control results

The **ELiTe InGenius software** interprets the PCR results for the target of the Positive Control and Negative Control reactions with the **ELiTe_PC** and **ELiTe_NC** Assay Protocols parameters. The resulting Ct values are converted to concentration and used to verify the system (reagents lot and instrument).

The Positive Control and Negative Control results, specific for the PCR reagent lot, are recorded in the database (Controls). They can be viewed and approved by "Administrator" or "Analyst" users following the GUI instructions.

The Positive Control and Negative Control results expire **after 15 days**.

The **ELiTe InGenius software** processes the Positive Control and Negative Control results and generates Control Charts. Four approved Positive Control and Negative Control results are used to set up the initial Control Chart. For subsequent controls, the results are analyzed by the software to ensure the system performances are within the acceptance criteria, shown in the Control Chart plots. Refer to the instrument manual for more details.

NOTE

If the Positive Control or Negative Control result does not meet the acceptance criteria, the "Failed" message is shown on the "Controls" screen. In this case, the results cannot be approved, and the Positive Control or Negative Control runs must be repeated.

NOTE

If the Positive Control or Negative Control result is not valid and samples were included in the same run, the samples can be approved but their results are not validated. In this case, the failed Control(s) and samples must all be repeated.

9.3.3 C. Validation of Sample results

The **ELiTe InGenius software** interprets the PCR results for the target (Channel **PJ**) and the Internal Control (Channel **IC**) with the **PJ ELiTe_BAL_200_100** and **PJ ELiTe_SP_200_100** Assay Protocol parameters. The resulting target Ct values are converted to concentration.

Results are shown in "Results Display" screen.

The sample results can be approved when the three conditions in the table below are true.

| 1) Calibration Curve | Status |
|----------------------|----------|
| PJ Q-PCR Standard | APPROVED |
| 2) Positive Control | Status |
| PJ Positive Control | APPROVED |

| 3) Negative Control | Status |
|---------------------|----------|
| PJ Negative Control | APPROVED |

The sample results are automatically interpreted by the **ELiTe InGenius software** using Assay Protocol parameters.

The possible result messages are listed in the table below.

For each sample the system reports a combination of the following messages specifying if the pathogen DNAs are either detected or not detected.

| Result of sample run | Interpretation |
|---|--|
| PJ: DNA Detected, quantity equal to XXX copies / mL | PJ DNA was detected in the sample within the assay measurement range, its concentration is shown. |
| PJ: DNA Detected, quantity below LLoQ copies / mL | PJ DNA was detected in the sample, its concentration is below the assay -Lower Limit of Quantification |
| PJ: DNA Detected, quantity beyond ULoQ copies / mL | PJ DNA was detected in the sample, its concentration is above the assay Upper Limit of Quantification |
| PJ: DNA Not detected or below LoD copies / mL | PJ DNA was not detected in the sample. The sample is negative for the target DNA, or its concentration is below the assay Limit of Detection . |
| Invalid - Retest Sample | Not valid assay result caused by Internal Control failure (due to e.g. incorrect extraction, inhibitors carry-over). The test should be repeated. |

Samples reported as “Invalid-Retest Sample”: in this case, the Internal Control DNA was not efficiently detected, which could be due to problems in sample collection, extraction or PCR steps (e. g. incorrect sampling, degradation or loss of DNA during the extraction or inhibitors in the eluate), which may cause incorrect results.

If sufficient eluate volume remains, the eluate can be retested (as is or diluted) by an amplification run in “PCR Only” mode. If the second result is invalid, the sample must be retested starting from extraction of a new sample using “Extract + PCR” mode (see [14 “Troubleshooting” page 29](#))

Samples reported as “PJ: DNA Not detected or below “LoD” copies/mL” are suitable for analysis but PJ DNA was not detected. In this case the sample may be either negative for PJ DNA or the PJ DNA is present at a concentration below the Limit of Detection of the assay (see [11 “Performance Characteristics” page 19](#)).

PJ DNA positive samples at a concentration below the Limit of Detection (and Lower Limit of Quantification) of the assay, if detected, are reported as “PJ: DNA Detected, quantity below “LLoQ” copies/mL” (see [11 “Performance Characteristics” page 19](#)).

PJ DNA positive samples within the Linear Measuring Range are detected and are reported as “PJ : DNA Detected, quantity equal to “XXX” copies / mL”. (see [11 “Performance Characteristics” page 19](#))

PJ DNA positive samples that are above the Upper Limit of Quantification are reported as “PJ: DNA Detected, quantity beyond “ULoQ” copies/mL” and they are not suitable for quantification. If needed the sample may be diluted before extraction or PCR and retested to yield results within the Linear Measuring Range of the assay (see [11 “Performance Characteristics” page 19](#)).

| NOTE |
|---|
| The results obtained with this assay must be interpreted in combination with all relevant clinical observation and laboratory outcomes. |

The sample results are stored in the database and, if valid, can be approved (Results Display) by “Administrator” or “Analyst” users, following the GUI instruction. From the “Results Display” window it is possible to print and save the Sample run results as “Sample Report” and “Track Report”.

9.3.4 D. Sample result reporting

The sample results are stored in the database and reports can be exported as “Sample Report” and “Track Report”.

The “Sample Report” shows the results details by selected sample (SID).

The “Track Report” shows the results details by selected Track.

The “Sample Report” and “Track Report” can be printed and signed by authorized personnel.

10 ELITE BeGenius PROCEDURE

The procedure to use the **Pneumocystis ELITE MGB Kit** with the **ELITE BeGenius** consists of three steps:

Table 7

| | | |
|--------|--------------------------------------|--|
| STEP 1 | Verification of the system readiness | |
| STEP 2 | Session setup | A) Sample run (Extract + PCR) |
| | | B) Eluted sample run (PCR Only) |
| | | C) Calibration run (PCR Only) |
| | | D) Positive Control and Negative Control run (PCR Only) |
| STEP 3 | Review and approval of results | A) Validation of Calibration curve |
| | | B) Validation of Positive Control and Negative Control results |
| | | C) Validation of sample results |
| | | D) Sample result reporting |

10.1 STEP 1 - Verification of the system readiness

Before starting the session:

- switch on the **ELITE BeGenius** and login in “**CLOSED**” mode,
- in the “Calibrations” menu on the Home page, verify the Calibrators (**Q - PCR Standard**) are approved and valid (Status) for the **PCR Mix** lot to be used. If no valid Calibrators are available for the **PCR Mix** lot, perform calibration as described in the following sections,
- in the “Controls” menu on the Home page, verify the PCR Controls (**Positive Control, Negative Control**) are approved and valid (Status) for the **PCR Mix** lot to be used. If no valid PCR Controls are available for the **PCR Mix** lot, run the PCR Controls as described in the following sections,

choose the type of run, following the instructions on the Graphical User Interface (GUI) for the session setup and using the Assay Protocols provided by EG SpA (see “Specimens and Controls”).

If the Assay Protocol of interest is not loaded in the system, contact your local ELITechGroup Customer Service.

10.2 STEP 2 – Session Setup

The **Pneumocystis ELITE MGB Kit** can be used on the **ELITE BeGenius** to perform:

- Sample run (Extract + PCR),
- Eluted sample run (PCR Only),
- Calibration run (PCR Only),
- Positive Control and Negative Control run (PCR Only).

All the required parameters are included in the Assay Protocol available on the instrument and are loaded automatically when the Assay Protocol is selected.

NOTE

The **ELiTe BeGenius** can be connected to the “Laboratory Information System” (LIS) which enables downloading the session information. Refer to the instrument manual for more details.

Before to setup a run:

Thaw the needed **PCR Mix** tubes at room temperature for 30 minutes. Each tube is sufficient for **12 tests** in optimized conditions (2 or more tests per session). Mix gently then spin down the contents for 5 seconds and keep on ice or cool block.

NOTE

Protect the **PCR Mix** from light while thawing because this reagent is photosensitive.

To set up one of the four types of run follow the steps below while referring to the GUI:

| | A. Sample run (Extract + PCR) | B. Eluted sample run (PCR Only) |
|----|--|---|
| 1 | Identify samples and, if needed, thaw at room temperature, mix gently, spin down the contents for 5 seconds and keep on ice or cool block. Thaw the needed CPE tubes at room temperature for 30 minutes. Mix gently, spin down the contents for 5 seconds and keep on ice or cool block. Each tube is sufficient for 12 extractions. | If needed, thaw the Elution tube containing the extracted nucleic acids at room temperature. Mix gently, then spin down the contents for 5 seconds and keep on ice or cool block. |
| 2 | Select “Perform Run” from the “Home” screen. | Select “Perform Run” from the “Home” screen. |
| 3 | Remove all the Racks from the “Cooler Unit” and place them on the preparation table. | Remove the “Racks” from “Lane 1, 2 and 3” (L1, L2, L3) of the “Cooler Unit” and place them on the preparation table. |
| 4 | Select the “Run mode”: “Extract + PCR”. | Select the “Run mode”: “PCR Only”. |
| 5 | Load the samples into the “Sample Rack”. (Note: when secondary tubes “2 mL Tubes” are loaded, use the blue adaptors for the “Sample Rack”). | Load the samples into the “Elution Rack”. |
| 6 | Insert the “Sample Rack” into the “Cooler Unit” starting from the “Lane 5” (L5). If needed, insert the “Sample ID” (SID) for each “Position” used. (If secondary tubes are loaded, flag “2 mL Tube”. If secondary tubes are not barcoded, type manually the “Sample ID”). | Insert the “Elution Rack” into the “Cooler Unit” starting from “Lane 3” (L3) If needed, for each “Position” enter the “Sample ID”, the “Sample matrix”, the “Extraction kit” and the “Extracted eluate vol.” (eluate volume). |
| 7 | Click “Next” to continue. | Click “Next” to continue. |
| 8 | Ensure the “Extraction Input Volume” is 200 µL and the “Extracted Elute Volume” is 100 µL. | Ensure the “Extraction Input Volume” is 200 µL and the “Extracted Elute Volume” is 100 µL. |
| 9 | Select the Assay Protocol in the “Assay” column (see “Specimens and Controls”). | Select the Assay Protocol in the “Assay” column (see “Specimens and Controls”). |
| 10 | Click “Next” to continue. | Click “Next” to continue. |
| 11 | When more than 12 samples are processed, repeat the procedure from point 6. | When more than 12 samples are processed, repeat the procedure from point 6. |
| 12 | Load the “Elution tubes” into the “Elution Rack” (Elution tubes can be labelled with barcode to improve traceability). | Not applicable |
| 13 | Insert the “Elution Rack” into the “Cooler Unit” starting from “Lane 3” (L3). When more than 12 samples are processed, repeat using “Lane 2” (L2). | Not applicable |

| | A. Sample run (Extract + PCR) | B. Eluted sample run (PCR Only) |
|----|---|--|
| 14 | Not applicable | Not applicable |
| 15 | Click "Next" to continue. | Not applicable |
| 16 | Load CPE and the PCR Mix into the "Reagent/Elution Rack". | Load the PCR Mix into "Reagent/Elution Rack". |
| 17 | Insert the "Reagent/Elution Rack" into the "Cooler Unit" in "Lane 2" (L2) if available or in "Lane 1" (L1). If needed, for each PCR Mix and / or CPE enter the "S/N" (serial number), the "Lot No." (lot number), the "Exp. Date" (expiry date) and the "T/R" (number of reactions). | Insert the "Reagent/Elution Rack" into the "Cooler Unit" in "Lane 2" (L2) if available or in "Lane 1" (L1). If needed, for each PCR Mix enter the "S/N" (serial number), the "Lot No." (lot number), the "Exp. Date" (expiry date) and the "T/R" (number of reactions). |
| 18 | Click "Next" to continue | Click "Next" to continue. |
| 19 | Verify the tips in the "Tip Rack(s)" in the "Inventory Area" and replace Tip Rack(s) if necessary. | Verify the tips in the "Tip Rack (s)" in the "Inventory Area" and replace Tip Rack(s) if necessary. |
| 20 | Click "Next" to continue. | Click "Next" to continue. |
| 21 | Load the "PCR Rack" with "PCR Cassette" in the Inventory Area. | Load the "PCR Rack" with "PCR Cassette" in the Inventory Area. |
| 22 | Click "Next" to continue. | Click "Next" to continue. |
| 23 | Load the "Extraction Rack" with the "ELITE InGenius SP 200" extraction cartridges and required extraction consumables. | Not applicable |
| 24 | Close the instrument door. | Close the instrument door. |
| 25 | Press "Start". | Press "Start". |

| | C. Calibration run (PCR Only) | D. Positive Control and Negative Control run (PCR Only) |
|---|---|---|
| 1 | Thaw the needed Q-PCR Standard tubes (Cal1: Q-PCR Standard 10 ² , Cal2: Q-PCR Standard 10 ³ , Cal3: Q-PCR Standard 10 ⁴ , Cal4: Q-PCR Standard 10 ⁵) for 30 minutes at room temperature. Mix gently then spin down the contents for 5 seconds and keep on ice or cool block. | Thaw the Positive Control tubes at room temperature for 30 minutes. Mix gently then spin down the contents for 5 seconds and keep on ice or cool block. Prepare the Negative Control by transferring at least 50 µL of molecular biology grade water to an "Elution tube", provided with the ELITE InGenius SP 200 Consumable Set. |
| 2 | Select "Perform Run" from the "Home" screen. | Select "Perform Run" from the "Home" screen |
| 3 | Remove the "Racks" from "Lane 1, 2 and 3" (L1, L2, L3) from the "Cooler Unit" and place them on the preparation table. | Remove the "Racks" from "Lane 1, 2 and 3" (L1, L2, L3) from the "Cooler Unit" and place them on the preparation table. |
| 4 | Select the "Run mode: PCR Only". | Select the "Run mode": "PCR Only". |
| 5 | Load the Q-PCR Standard tubes into the "Elution Rack". | Load the Positive Control and Negative Control tubes into the "Elution Rack". |
| 6 | Insert the " Elution Rack " into the "Cooler Unit" starting from the "Lane 3" (L3). If needed, for each "Position" enter the "Reagent name" and the "S/N" (serial number), the "Lot No." (lot number), the "Exp. Date" (expiry date) and the "T/R" (number of reactions). | Insert the " Elution Rack " into the "Cooler Unit" starting from the "Lane 3" (L3). If needed, for each "Position" enter the "Reagent name" and the "S/N" (serial number), the "Lot No." (lot number), the "Exp. Date" (expiry date) and the "T/R" (number of reactions). |
| 7 | Click "Next" to continue. | Click "Next" to continue. |
| 8 | Ensure the "Extraction Input Volume" (200 µL) and the "Extracted Elute Volume" (100 µL). | Ensure the "Extraction Input Volume" (200 µL) and the "Extracted Elute Volume" (100 µL). |

| | C. Calibration run (PCR Only) | D. Positive Control and Negative Control run (PCR Only) |
|----|--|---|
| 9 | Select the Assay Protocol in the “Assay” column (see “Specimens and Controls”). | Select the Assay Protocol in the “Assay” column (see “Specimens and Controls”). |
| 10 | Click “Next” to continue. | Click “Next” to continue. |
| 11 | Load the PCR Mix into “Reagent/Elution Rack”. | Load the PCR Mix into “Reagent/Elution Rack”. |
| 12 | Insert the “ Reagent/Elution Rack ” into the “Cooler Unit” in “Lane 2” (L2) If needed, for each PCR Mix enter the “S/N” (serial number), the “Lot No.” (lot number), the “Exp. Date” (expiry date) and the “T/R” (number of reactions). | Insert the “Reagent/Elution Rack” into the “Cooler Unit” in “Lane 2” (L2). If needed, for each PCR Mix enter the “S/N” (serial number), the “Lot No.” (lot number), the “Exp. Date” (expiry date) and the “T/R” (number of reactions). |
| 13 | Click “Next” to continue. | Click “Next” to continue. |
| 14 | Verify the tips in the “Tip Racks” in the “Inventory Area” and replace Tip Rack(s) if necessary. | Verify the tips in the “Tip Racks” in the “Inventory Area” and replace Tip Rack(s) if necessary. |
| 15 | Click “Next” to continue. | Click “Next” to continue. |
| 16 | Load the “ PCR Rack ” with “ PCR Cassette ” in the Inventory Area. | Load the “ PCR Rack ” with “ PCR Cassette ” in the Inventory Area. |
| 17 | Click “Next” to continue. | Click “Next” to continue. |
| 18 | Close the instrument door. | Close the instrument door. |
| 19 | Press “Start”. | Press “Start”. |

When the session is finished, the **ELITE BeGenius** allows users to view, approve, store the results, print and save the report.

NOTE

At the end of the run, the remaining Extracted Sample in the **Elution tube** must be removed from the instrument, capped, identified, and stored at $-20 \pm 10^{\circ}\text{C}$ for no longer than one month. Avoid the spilling of the Extracted Sample.

NOTE

At the end of the run the **PCR Mix** can be removed from the instrument, capped and stored at -20°C or below or can be kept on board in the refrigerated block for up to 7 hours (2 sessions of 3 hours each and the time needed to start a third session); mix gently and spin down the content for 5 seconds before starting the next session.

NOTE

At the end of the run, the remaining **Q - PCR Standard** can be removed from the instrument, capped and stored at -20°C or below. Avoid spilling the Q - PCR Standard.

NOTE

The **Q- PCR Standard** can be used for 4 separate sessions of 2 hours each.

NOTE

At the end of the run, the remaining **Positive Control** can be removed from the instrument, capped and stored at -20°C or below. Avoid the spilling of the **Positive Control**. The remaining **Negative Control** must be discarded.

NOTE

The **Positive Control** can be used for 4 separate sessions of 3 hours each.

NOTE

At the end of the run, the **PCR Cassette** and the other consumables must be disposed of following all governmental and environmental regulations. Avoid spilling the reaction products.

10.3 STEP 3 -Review and approval of results

The **ELiTe BeGenius** monitors target and internal control fluorescence signals for each reaction and automatically applies the Assay Protocol parameters to generate PCR curves which are then interpreted into results.

At the end of the run, the “Results Display” screen is automatically shown. In this screen the results and the run information are shown. From this screen results can be approved, and reports printed or saved (“Sample Report” or “Track Report”). Refer to the instrument manual for more details.

NOTE

The **ELiTe BeGenius** can be connected to the “Laboratory Information System” (LIS) which enables uploading the session results to the laboratory data center. Refer to the instrument manual for more details.

The **ELiTe BeGenius** generates the results with the **Pneumocystis ELiTe MGB Kit** through the following procedure:

- A. Validation of Calibration curve,
- B. Validation of Positive Control and Negative Control results,
- C. Validation of sample results,
- D. Sample result reporting.

NOTE

Please, refer to the same paragraph of the **ELiTe InGenius Procedure** for the details.

11 PERFORMANCE CHARACTERISTICS**11.1 Analytical sensitivity: Limit of Detection (LoD)**

The analytical sensitivity of the Pneumocystis ELiTe MGB Kit, as Limit of Detection (LoD) was defined in association with BAL/BA samples and **ELiTe InGenius** system.

The LoD was calculated by testing a panel of negative BAL/BA samples spiked with *Pneumocystis jirovecii* (PJ) certified reference material at known titre (Qnostics). The LoD was obtained by Probit regression analysis of the data as the concentration corresponding to 95% probability of a positive call.

The final results are reported in the following table.

Table 8 Limit of Detection in BAL/BA with ELiTe InGenius

| LoD | 95% confidence interval | |
|----------------|-------------------------|-----------------|
| | lower bound | upper bound |
| 97 copies / mL | 60 copies / mL | 275 copies / mL |

The calculated LoD value was verified for each matrix by testing on **ELiTe InGenius** and on **ELiTe BeGenius** a pool of matrix spiked by Pneumocystis reference material “at the claimed concentration.

The results obtained confirmed the claimed concentration for the target of Pneumocystis ELITE MGB Kit on both ELITE InGenius and ELITE BeGenius for each matrix.

11.2 Linear measuring range

The linear measuring range of the assay was determined in association with each matrix on **ELITE InGenius** and **ELITE BeGenius** using a panel of dilution of PJ target DNA (spiked with a quantified plasmid DNA containing the amplicon of the mtLSU) in PJ DNA - negative matrix.

The results, for each matrix, are reported in the following paragraphs.

Sputum:

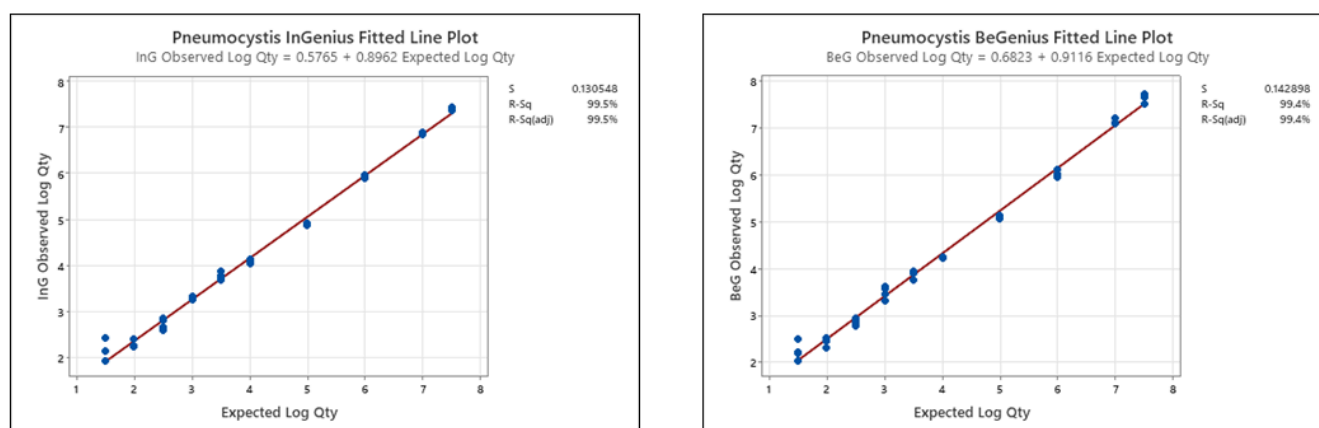


Fig. 1

The results obtained by ELITE InGenius and ELITE BeGenius were analyzed by orthogonal and linear regression in order to calculate the correlation between the methods.

The results are summed up in the following figure.

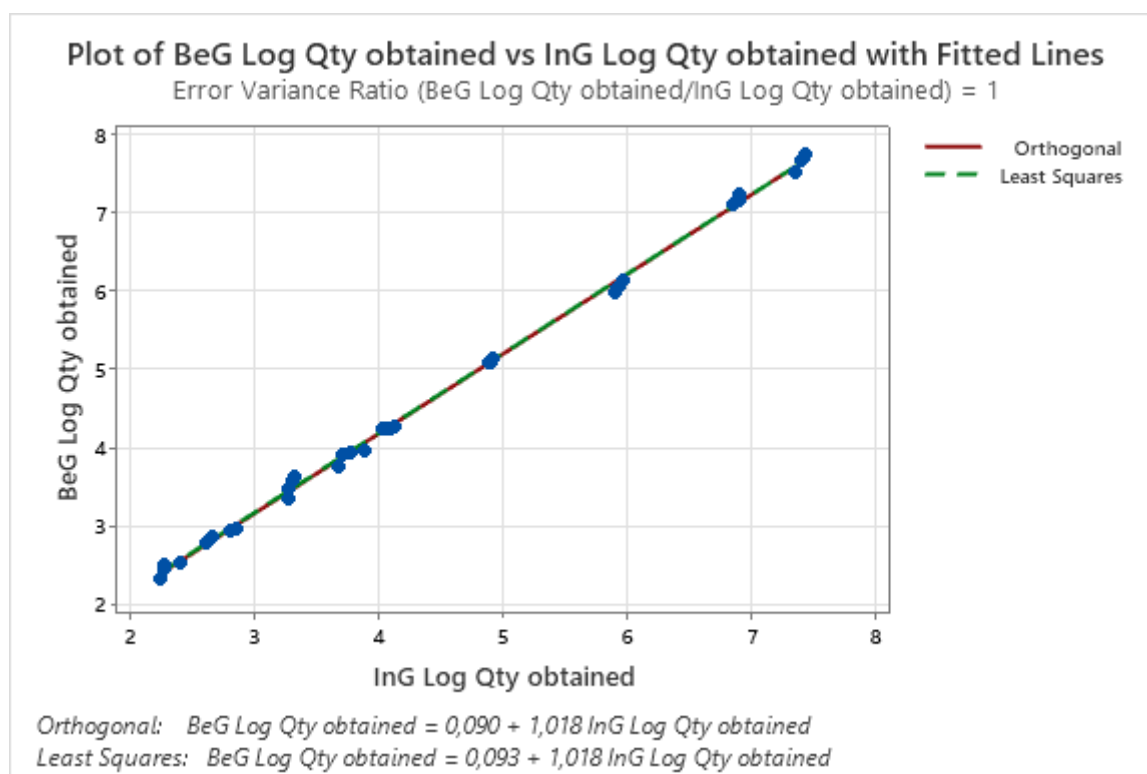


Fig. 2

The Orthogonal Regression analysis generated an intercept equal to 0.090 (95% CI: 0.0291 0.1509) and a slope equal to 1.018 (95% CI: 1.0059, 1.0308).

BAL/BA:

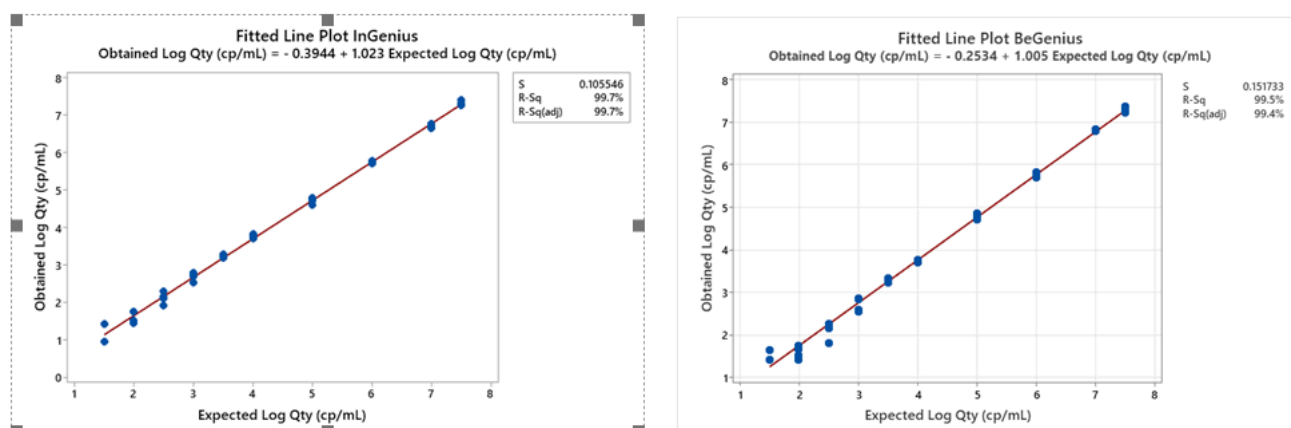


Fig. 3

The results obtained by ELITe InGenius and ELITe BeGenius were analyzed by orthogonal and linear regression in order to calculate the correlation between the methods.

The results are summed up in the following figure.

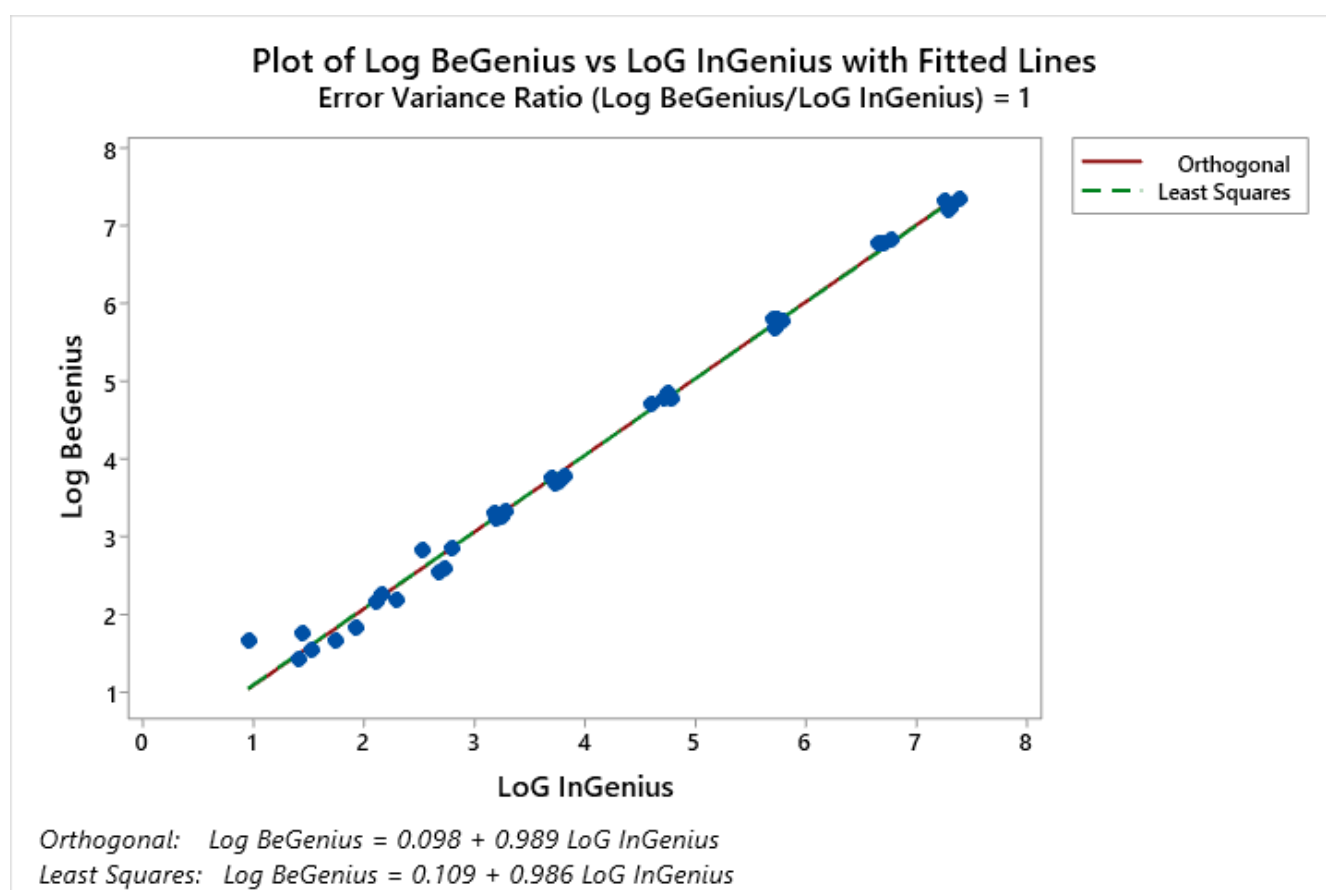


Fig. 4

The Orthogonal Regression analysis generated an intercept equal to 0.098 (95% CI -0.0134, 0.2099) and a slope equal to 0.989 (95% CI: 0.9647, 1.0137).

The final results are summarized in the following table.

Table 9 Linear measuring range for Sputum and BAL/BA and ELiTe InGenius and ELiTe BeGenius

| Lower Limit | Upper limit |
|----------------|------------------------|
| 97 copies / mL | 31,622,777 copies / mL |

11.3 Inclusivity: Efficiency of detection and quantification efficiency on different genotypes

The efficiency of detection for different genotypes was evaluated by comparison with sequences available in the EBI ENA nucleotide database.

The analysis of the regions chosen for the hybridization of the primers and of the fluorescent probes in the alignment of the sequences available in the database for the region of the mtLSU gene of PJ isolates showed their conservation and absence of significant mutations.

The inclusivity of the assay, as detection efficiency of different isolates of *Pneumocystis* was evaluated by testing the QCMD 2018 *Pneumocystis jirovecii pneumonia* EQA Panel (Qnostics) composed with 6 different clinical isolates.

The final results are reported in the following table.

Table 10

| Reference material | | | | Pneumocystis ELiTe MGB Kit results | |
|--------------------|--------------------------------|---------------------|---------------|------------------------------------|--------------|
| Samples | Sample description | Detection Frequency | Sample Status | PJ Ct | PJ copies/mL |
| PCPDNA18S-01 | <i>P. jirovecii</i> (Clinical) | Detected | Core | 33.60 | 1,771 |
| | | | | 33.69 | 2,704 |
| PCPDNA18S-02 | <i>P. jirovecii</i> (j888023) | Frequently Detected | Core | 32.19 | 4,805 |
| | | | | 32.56 | 5,706 |
| PCPDNA18S-03 | <i>P. jirovecii</i> (j888023) | Frequently Detected | Core | 30.11 | 20,904 |
| | | | | 30.28 | 25,449 |
| PCPDNA18S-04 | <i>P. jirovecii</i> (Type A1) | Frequently Detected | Core | 33.83 | 1,507 |
| | | | | 33.57 | 2,933 |
| PCPDNA18S-05 | <i>P. jirovecii</i> (Clinical) | Detected | Core | 34.06 | 1,280 |
| | | | | 34.81 | 1,298 |
| PCPDNA18S-06 | <i>P. jirovecii</i> (g885652) | Frequently Detected | Core | 35.60 | 432 |
| | | | | 37.69 | 196 |
| PCPDNA18S-07 | <i>P. jirovecii</i> (j888023) | Detected | Educational | 37.55 | 109 |
| | | | | 37.14 | 280 |
| PCPDNA18S-08 | <i>P. jirovecii</i> Negative | Negative | Core | Undet | Undet |
| | | | | Undet | Undet |
| PCPDNA18S-09 | <i>P. jirovecii</i> (Type A1) | Detected | Educational | 37.11 | 149 |
| | | | | 36.91 | 326 |

Table 10 (continued)

| | | | | | |
|--------------|----------------------------------|------------------------|------|-------|-------|
| PCPDNA18S-10 | <i>P. jirovecii</i> (j888023) | Frequently Detected | Core | 31.75 | 6,547 |
| | | | | 33.27 | 3,565 |

All of the 9 positive samples were correctly detected by the Pneumocystis ELITE MGB Kit in association with the ELITE InGenius instrument.

The inclusivity of the assay, as detection efficiency of different isolates of Pneumocystis was also evaluated by analyzing positive clinical samples of BAL/BA and Sputum tested by a reference CE IVD method.

The results are summarized in the following table.

Table 11

| Samples | N | positive | negative | invalid |
|--------------------|----|----------|----------|---------|
| PJ positive BAL/BA | 15 | 15 | 0 | 0 |
| PJ positive Sputum | 8 | 8 | 0 | 0 |

All of the 23 positive clinical samples were found positive by the Pneumocystis ELITE MGB Kit in association with the ELITE InGenius instrument.

11.4 Potential interfering markers: Cross-reactivity

The potential cross-reactivity with other unintended organisms of the product Pneumocystis ELITE MGB Kit was evaluated by *in silico* analysis of sequences available in the EBI ENA nucleotide database. The analysis showed no significant homology with other unintended organisms (other organisms, including other fungi).

The absence of cross-reactivity with potential interfering organisms that can be found in clinical samples of BAL/BA was also verified through the analysis of a panel of certified reference materials (ATCC) at high titer.

The final results are reported in the following table.

Table 12 Potential cross-reactive markers

| Organism | Strain | Outcome |
|-----------------------------------|----------------|---------------------|
| <i>Aspergillus fumigatus</i> | 118 | No cross-reactivity |
| <i>Candida albicans</i> | 3147 | No cross-reactivity |
| <i>Staphylococcus aureus</i> | Rosenbach | No cross-reactivity |
| <i>Escherichia coli</i> | H10407 | No cross-reactivity |
| <i>Bordetella pertussis</i> | - | No cross-reactivity |
| <i>Bordetella parapertussis</i> | - | No cross-reactivity |
| <i>Haemophilus influenzae</i> | Rd | No cross-reactivity |
| <i>Streptococcus pneumoniae</i> | R6 | No cross-reactivity |
| <i>Legionella pneumophila</i> | Philadelphia-1 | No cross-reactivity |
| <i>Mycoplasma pneumoniae</i> | FH | No cross-reactivity |
| <i>Chlamydomonas pneumoniae</i> | AR-39 | No cross-reactivity |
| <i>Mycobacterium tuberculosis</i> | - | No cross-reactivity |

Table 12 Potential cross-reactive markers (continued)

| | | |
|-------------------|------------------|---------------------|
| CMV | AD-169 | No cross-reactivity |
| Enterovirus | Enterovirus 71 | No cross-reactivity |
| Adenovirus | Adenoid 6 | No cross-reactivity |
| Influenza A virus | A/PR/8/34 | No cross-reactivity |
| Influenza B virus | B/Florida/4/2006 | No cross-reactivity |
| RSV | A2 | No cross-reactivity |

All the potential cross-reacting organisms gave a negative result for the target when tested by Pneumocystis ELITE MGB Kit.

11.5 Potential interfering markers: Inhibition

The absence of interference on target amplification by other organisms that can be found in BAL/BA clinical samples was verified by testing the panel of certified materials. Genomic DNA or RNA from different potentially cross-reactive organisms (ATCC) were spiked with plasmid DNA containing the amplicon of the mtLSU at low concentration (about 10 copies / reaction).

The final results are reported in the following table.

Table 13 Potential interfering markers

| Organism | Strain | Outcome |
|-----------------------------------|------------------|-----------------|
| <i>Aspergillus fumigatus</i> | 118 | No interference |
| <i>Candida albicans</i> | 3147 | No interference |
| <i>Staphylococcus aureus</i> | Rosenbach | No interference |
| <i>Escherichia coli</i> | H10407 | No interference |
| <i>Bordetella pertussis</i> | - | No interference |
| <i>Bordetella parapertussis</i> | - | No interference |
| <i>Haemophilus influenzae</i> | Rd | No interference |
| <i>Streptococcus pneumoniae</i> | R6 | No interference |
| <i>Legionella pneumophila</i> | Philadelphia-1 | No interference |
| <i>Mycoplasma pneumoniae</i> | FH | No interference |
| <i>Chlamydia pneumoniae</i> | AR-39 | No interference |
| <i>Mycobacterium tuberculosis</i> | - | No interference |
| CMV | AD-169 | No interference |
| Enterovirus | Enterovirus 71 | No interference |
| Adenovirus | Adenoid 6 | No interference |
| Influenza A virus | A/PR/8/34 | No interference |
| Influenza B virus | B/Florida/4/2006 | No interference |
| RSV | A2 | No interference |

All the potential interfering organisms did not impair the amplification of the target when tested by Pneumocystis ELITE MGB Kit.

11.6 Potential interfering substances

The possible effect of interfering substances was evaluated with samples containing endogenous substances: whole blood, mucin, antibiotics, steroids, antihistamines.

BAL/BA samples were spiked with reference material of PJ (Qnostics) at a concentration of 3 x LoD and with the potential interfering substances.

The Ct values of the target and the Internal Control (reference and test samples) were used to calculate the percentage Coefficient of Variability (%CV) in order to evaluate the possible interference.

The results are reported in the following table.

Table 14 Potential interfering substances

| Sample | Pos. / Rep. | %CV PJ Ct | %CV IC Ct |
|--|-------------|-----------|-----------|
| 5% Whole Blood | 3 / 3 | 1.41 | 0.94 |
| 1% Mucin | 3 / 3 | 0.29 | 0.91 |
| 6 µg / mL Ambroxol Hydrochloride | 3 / 3 | 0.60 | 0.53 |
| 160 µg / mL Sulfamethoxazole + 32 µg / mL Trimethoprim | 3 / 3 | 1.35 | 0.54 |
| 200 µg / mL Ampicillin | 3 / 3 | 1.72 | 0.60 |
| 40 ng / mL Beclomethasone | 3 / 3 | 1.28 | 0.71 |
| 4 µg / mL Ebastine | 3 / 3 | 0.64 | 0.81 |

All the samples resulted positive for the target of interest as expected. The percentage %CV of Ct values were lower than 2%. None of the tested substances at the tested concentrations were found to interfere with the target detection by Pneumocystis ELITE MGB Kit.

11.7 Repeatability

The Repeatability of the assay was evaluated on **ELITE InGenius** and **ELITE BeGenius** instruments by analysis of a panel of BAL/BA samples, including one negative sample and two samples spiked with PJ DNA.

Results are shown in the tables below.

Table 15 Intra - Session Repeatability on ELITE InGenius

| Sample | N | PJ | | | |
|----------|---|---------|------|------|-------------|
| | | Mean Ct | SD | % CV | % Agreement |
| Negative | 8 | NA | NA | NA | 100% |
| 3X LoD | 8 | 37.08 | 0.42 | 1.13 | 100% |
| 10X LoD | 8 | 35.04 | 0.48 | 1.37 | 100% |

Table 16 Inter - Session Repeatability on ELiTe InGenius

| Sample | N | PJ | | | |
|----------|----|---------|------|------|-------------|
| | | Mean Ct | SD | % CV | % Agreement |
| Negative | 16 | NA | NA | NA | 100% |
| 3X LoD | 16 | 37.23 | 0.72 | 1.93 | 100% |
| 10X LoD | 16 | 35.28 | 0.46 | 1.30 | 100% |

Table 17 Intra - Session Repeatability on ELiTe BeGenius

| Sample | N | PJ | | | |
|----------|---|---------|------|------|-------------|
| | | Mean Ct | SD | % CV | % Agreement |
| Negative | 8 | NA | NA | NA | 100% |
| 3X LoD | 8 | 37.45 | 0.32 | 0.86 | 100% |
| 10X LoD | 8 | 35.37 | 0.44 | 1.25 | 100% |

Table 18 Inter - Session Repeatability on ELiTe BeGenius

| Sample | N | PJ | | | |
|----------|----|---------|------|------|-------------|
| | | Mean Ct | SD | % CV | % Agreement |
| Negative | 16 | NA | NA | NA | 100% |
| 3X LoD | 16 | 37.98 | 1.05 | 2.76 | 100% |
| 10X LoD | 16 | 35.79 | 0.63 | 1.76 | 100% |

In the Repeatability test, the Pneumocystis ELiTe MGB Kit detected all the samples as expected and showed a maximum variability of target Ct values as %CV lower than 5%.

11.8 Reproducibility

The Reproducibility of the assay was evaluated on **ELiTe InGenius** and **ELiTe BeGenius** instruments by analysis of a panel of BAL/BA samples, including one negative sample and two samples spiked with PJ DNA.

A summary of Inter-Instrument Reproducibility is shown in the tables.

Table 19 Inter – Instrument Reproducibility on ELiTe InGenius

| Sample | N | PJ | | | |
|----------|---|---------|------|------|-------------|
| | | Mean Ct | SD | % CV | % Agreement |
| Negative | 8 | NA | NA | NA | 100% |
| 3X LoD | 8 | 36.63 | 0.63 | 1.72 | 100% |
| 10X LoD | 8 | 34.42 | 0.48 | 1.39 | 100% |

Table 20 Inter – Instrument Reproducibility on ELiTe BeGenius

| Sample | N | PJ | | | |
|----------|---|---------|------|------|-------------|
| | | Mean Ct | SD | % CV | % Agreement |
| Negative | 8 | NA | NA | NA | 100% |
| 3X LoD | 8 | 37.26 | 0.41 | 1.11 | 100% |
| 10X LoD | 8 | 34.92 | 0.40 | 1.16 | 100% |

Table 21 Inter – Batch Reproducibility on ELiTe InGenius

| Sample | N | PJ | | | |
|----------|---|---------|------|------|-------------|
| | | Mean Ct | SD | % CV | % Agreement |
| Negative | 8 | NA | NA | NA | 100% |
| 3X LoD | 8 | 36.99 | 0.48 | 1.29 | 100% |
| 10X LoD | 8 | 35.03 | 0.44 | 1.25 | 100% |

Table 22 Inter – Batch Reproducibility on ELiTe BeGenius

| Sample | N | PJ | | | |
|----------|---|---------|------|------|-------------|
| | | Mean Ct | SD | % CV | % Agreement |
| Negative | 8 | NA | NA | NA | 100% |
| 3X LoD | 8 | 37.59 | 0.38 | 1.02 | 100% |
| 10X LoD | 8 | 35.70 | 0.59 | 1.65 | 100% |

In the Reproducibility test, the Pneumocystis ELiTe MGB Kit detected all the samples as expected and showed a maximum variability of target Ct values as %CV lower than 5%.

11.9 Diagnostic specificity: confirmation of negative samples

The Diagnostic specificity of the assay, as confirmation of negative clinical samples was evaluated in association with **ELiTe InGenius** by analyzing BAL/BA and Sputum samples that were negative for PJ DNA (tested with a CE IVD real time amplification product). As **ELiTe BeGenius** has equivalent analytical performances to **ELiTe InGenius**, the diagnostic performances of the assay performed on the two instruments are also considered equivalent. Therefore, the Diagnostic specificity of the assay obtained in association with **ELiTe InGenius** is also applicable to **ELiTe BeGenius**.

The results are summarized in the following table.

| Samples | N | Positive | Negative | %Diagnostic Specificity |
|--------------------|----|----------|----------|-------------------------|
| PJ Negative BAL/BA | 58 | 2 | 56 | 96.6 |
| PJ Negative Sputum | 35 | 1 | 34 | 97.10 |

All samples were valid for analysis.

The Internal Control Ct (IC Ct) cut-off value is set at 34 with **ELiTe InGenius** and at 35 with **ELiTe BeGenius** for BAL/BA samples.

The Internal Control Ct (IC Ct) cut-off value is set at 37 for sputum samples when tested with **ELiTe InGenius** and **ELiTe BeGenius**.

11.10 Diagnostic sensitivity: confirmation of positive samples

The Diagnostic sensitivity of the assay, as confirmation of positive clinical samples was evaluated in association with ELiTe InGenius by analyzing clinical BAL/BA and Sputum samples positive for PJ DNA (tested with a CE IVD real time PCR kit) or spiked with PJ reference material. As ELiTe BeGenius has equivalent analytical performances to ELiTe InGenius, the diagnostic performances of the assay performed on the two instruments are also considered equivalent. Therefore, the Diagnostic sensitivity of the assay obtained in association with ELiTe InGenius is also applicable to ELiTe BeGenius.

The results are summarized in the following table

| Samples | N | Positive | Negative | %Diagnostic Sensitivity |
|-----------------------------|-----------|-----------|----------|-------------------------|
| PJ positive BAL/BA | 13 | 13 | 0 | 100 |
| PJ spiked BAL/BA | 30 | 30 | 0 | 100 |
| Total BAL/BA samples | 43 | 43 | 0 | 100 |
| PJ positive Sputum | 8 | 8 | 0 | 100 |
| PJ spiked Sputum | 30 | 24 | 6 | 80 |
| Total Sputum samples | 38 | 32 | 6 | 84.2 |

NOTE

The complete data and results of the tests carried out to evaluate the product performance characteristics with matrices and instrument are recorded in the Product Technical File "PJ ELiTe MGB Kit", FTP150ING.

12 REFERENCES

- C. Valero et al. (2016) *Front. Microbiol.* 7:1413
M. Maillet et al. (2014) *Eur J Clin Microbiol Infect Dis.* 33(3):331-6
E. A. Lukhtanov et al. (2007) *Nucleic Acids Res.* 35: e30
K. Linnet et al. (2004) *Clin. Chem.* 50: 732 - 740

13 PROCEDURE LIMITATIONS

Use this product only with the following clinical sample: BAL/BA, liquefied Sputum.

Currently there are no data available concerning product performance with other clinical samples.

Do not use this product with samples containing too much mucin: samples with high viscosity inhibit the amplification reaction of nucleic acids and can cause invalid results.

There are no data available concerning product performance with DNA extracted from the following clinical samples: respiratory swabs.

The results obtained with this product depend on proper identification, collection, transport storage and processing of the samples. To avoid incorrect results, it is therefore necessary to take care during these steps and to carefully follow the instructions for use provided with the product.

Owing to its high analytical sensitivity, the Real Time PCR method used in this product is sensitive to contamination from positive clinical samples, positive controls and PCR products. Cross-contamination cause false positive results. The product format is designed to limit cross-contamination. However, cross-contamination can only be avoided by good laboratory practices and following these instructions for use.

This product must be handled by qualified personnel trained in the processing of potentially infective biological samples and chemical preparations classified as dangerous to prevent accidents with potentially serious consequences for the user and other persons.

This product requires the use of personal protective equipment and areas that are suitable for the processing of potentially infective biological samples and chemical preparations classified as dangerous to prevent accidents with potentially serious consequences for the user and other persons.

This product requires the use of personal protective equipment and instruments dedicated to work session setup to avoid false positive results.

To avoid incorrect results, this product must be handled by professional personnel, qualified and trained in molecular biology techniques such as extraction, PCR and detection of nucleic acids.

Due to inherent differences between technologies, it is recommended that users perform method correlation studies to estimate technology differences prior to switching to a new technology.

A negative result obtained with this product indicates that the target DNA is not detected in the DNA extracted from the sample; however, it cannot be excluded that the target DNA has a lower titer than the product detection limit (see Performance Characteristics). In this case the result could be a false negative.

Results obtained with this product may sometimes be invalid due to failure of internal control. In this case the sample shall be retested, starting from extraction, which can lead to a delay in obtaining final results.

Possible polymorphisms, insertions or deletions within the region of the DNA targeted by the product primers and probes may impair detection and quantification of target DNA.

As with any other diagnostic medical device, the results obtained with this product must be interpreted in combination with all relevant clinical observations and laboratory results.

As with any other diagnostic medical device, there is a residual risk of obtaining invalid, or erroneous results with this product. This residual risk cannot be eliminated or further reduced. In some cases, this residual risk could contribute to wrong decisions with potentially dangerous effects for the patient.

However, this residual risk associated to the intended use of the product has been weighed against the potential benefits to the patient and it has been assessed acceptable.

14 TROUBLESHOOTING

Table 23

| Invalid Q-PCR Standard reaction, Standard curve or Positive Control reaction | |
|--|--|
| Possible Causes | Solutions |
| Instrument setting error. | Check the position of PCR Mix, Q-PCR Standards and Positive Control. Check the volumes of PCR Mix, Q-PCR Standards and Positive Control. |
| PCR Mix degradation. | Do not use the PCR Mix for more than 7 independent sessions (3 hours each in the Inventory Area Cool Block or in the Cooler Unit). Do not use the PCR Mix for more than three consecutive sessions (7 hours in the Inventory Area Cool Block or in the Cooler Unit). Do not leave the PCR Mix at room temperature for more than 30 minutes. Use a new aliquot of PCR Mix. |
| Q-PCR Standards or Positive Control degradation. | Do not use the Q-PCR Standard for more than 4 independent sessions (2 hours each in the Extraction Area or in the Cooler Unit). Do not use the Positive Control for more than 4 independent sessions (3 hours each in the Extraction Area or in the Cooler Unit). Use new aliquots of Q-PCR Standards or Positive Control. |
| Instrument error. | Contact ELiTechGroup Technical Service. |

Table 24

| Invalid Negative Control reaction | |
|---|---|
| Possible Causes | Solutions |
| Instrument setting error. | Check the position of PCR Mix and Negative Control. Check the volumes of PCR Mix and Negative Control. |
| Contamination of the PCR Mix. | Use a new aliquot of PCR Mix. |
| Contamination of the Negative Control. | Do not use the Negative Control for more than one session. Use a new aliquot of molecular biology grade water. |
| Contamination of the extraction area, of Racks or of Inventory Block. | Clean surfaces with aqueous detergents, wash lab coats, replace tubes and tips in use. |
| Instrument error. | Contact ELITechGroup Technical Service. |

Table 25

| Invalid Sample reaction | |
|---|---|
| Possible Causes | Solutions |
| Instrument setting error. | Check the position of PCR Mix, Internal Control, and sample. Check the volumes of PCR Mix, Internal Control, and sample. |
| PCR Mix degradation. | Do not use the PCR Mix for more than 7 independent sessions (3 hours each in the Inventory Area or in the Cooler Unit). Do not use the PCR Mix for more than three consecutive sessions (7 hours in the Inventory Area Cool Block or in the Cooler Unit). Do not leave the PCR Mix at room temperature for more than 30 minutes. Use a new aliquot of PCR Mix. |
| Internal Control template degradation. | Use a new aliquot of Internal Control. |
| Inhibition due to interfering substances in the sample. | Repeat the amplification with a 1:2 dilution in molecular biology grade water of eluted sample in a "PCR Only" session. Repeat the extraction with a 1:2 dilution in molecular biology grade water of the sample in an "Extract + PCR" session. |
| Instrument error. | Contact ELITechGroup Technical Service. |

Table 26

| Anomalous dissociation curve | |
|---|--|
| Possible causes | Solutions |
| Absence of a defined peak. Defined peak but T _m different from that of the other samples and that of the Standards or Positive Control. | Check for target Ct lower than 30. High quantity of amplification product at the end of the reaction may interfere with the melting curve analysis. Repeat the sample amplification to confirm the presence of target with a possible mutation. The target in the sample should be sequenced to confirm mutation. |

Table 27

| Error in Ct calculation | |
|---|--|
| Possible Causes | Solutions |
| Too high concentration of target in the sample with anomalous fluorescence signal | <p>If significant amplification is observed in PCR plot select the track related to the sample and manually approve the result as positive.</p> <p>If no amplification is observed in PCR plot select the track related to the sample and manually approve the result as negative or leave it as invalid.</p> <p>If a Ct value is required:</p> <ul style="list-style-type: none"> - repeat the amplification of eluted sample with a 1:10 dilution in molecular biology grade water in a "PCR Only" session or - repeat the extraction of the sample with a 1:10 dilution in molecular biology grade water in an "Extract + PCR" session. |

Table 28

| Abnormal high rate of positive results within the same session (reactions with similar late Ct values) | |
|--|---|
| Possible Causes | Solutions |
| Sample-to-sample contamination during preanalytical steps. | <p>Clean the micropipette with fresh 3% sodium hypochlorite solution (bleach) or DNA/RNA cleaner after pipetting each sample.</p> <p>Do not use Pasteur pipettes. The pipettes must be of the positive displacement type or used with aerosol filter tips.</p> <p>Introduce samples in the last positions of the instruments, as indicated by the GUI. Follow the loading sequence indicated by the software.</p> |
| Laboratory environmental contamination. | <p>Clean all surfaces in contact with the operator and samples (including the pipettes) with fresh 3% sodium hypochlorite solution (bleach) or DNA/RNA cleaner.</p> <p>Perform an U.V. decontamination cycle.</p> <p>Use a new tube of PCR Mix and / or CPE.</p> |

15 SYMBOLS



Catalogue Number.



Upper limit of temperature.



Batch code.



Use by (last day of month).



in vitro diagnostic medical device.



Fulfilling the requirements of the European Directive 98/79/EC for *in vitro* diagnostic medical device.



Unique Device Identification



Contains sufficient for "N" tests.



Caution, consult instructions for use.



Contents.



Keep away from sunlight.



Manufacturer.

16 NOTICE TO PURCHASER: LIMITED LICENSE

This product contains reagents manufactured by Thermo Fisher Scientific and are sold under licensing arrangements between ELiTechGroup S.p.A. and its Affiliates and Thermo Fisher Scientific. The purchase price of this product includes limited, nontransferable rights to use only this amount of the product solely for activities of the purchaser which are directly related to human diagnostics. For information on purchasing a license to this product for purposes other than those stated above, contact Licensing Department, Thermo Fisher Scientific. Email: outlicensing@thermofisher.com.

ELiTe MGB Kit® detection reagents are covered by one or more of U.S. Patent numbers 7319022, 7348146, 7381818, 7541454, 7671218, 7718374, 7723038, 7759126, 7767834, 8008522, 8067177, 8163910, 8389745, 8969003, 9056887, 9085800, 9169256, 9328384, 10677728, 10738346, 10890529, and EP patent numbers 1687609, 1781675, 1789587, 2689031, 2714939, 2736916, 2997161 as well as applications that are currently pending.

ELiTe InGenius® and ELiTe BeGenius® technologies are covered by patents and pending applications.

This limited license allows the person or entity to whom the product has been provided to use the product and data generated by the use of the product, solely for human diagnostics. Neither ELiTechGroup S.p.A. nor its licensors grant any other licenses, expressed or implied for any other purposes.

Appendix A Pneumocystis ELiTe MGB Kit used in association with Genius series® platforms



CAUTION

This document is a simplified version of the official instruction for use. Please refer to the complete document before use: www.elitechgroup.com

INTENDED USE

The product **Pneumocystis ELiTe MGB Kit** is an *in vitro* diagnostic medical device intended to be used by healthcare professionals as a quantitative nucleic acids Real-Time PCR assay for the **detection and the quantification of the DNA of *Pneumocystis jirovecii* (PJ)** in DNA samples extracted from clinical specimens.

The assay is validated in association with the **ELiTe InGenius®** and **ELiTe BeGenius®** instruments, automated and integrated systems for extraction, Real-Time PCR and results interpretation, using human Broncho-alveolar lavage (BAL) / bronchial aspirate (BA) and liquefied Sputum specimens

The product is intended for *in vitro* diagnostic use as an aid in diagnosis and monitoring of *Pneumocystis jirovecii* infection.

The results must be interpreted in combination with all relevant clinical observations and laboratory outcomes.


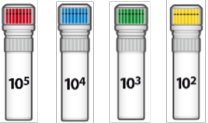

Amplified sequence

| Sequence | Gene | Fluorophore | Channel |
|------------------|---------------------------|-------------|---------|
| Target | mtLSU | FAM | PJ |
| Internal Control | artificial sequence (IC2) | AP525 | IC |

Validated matrices

- BAL/BA
- Liquefied Sputum samples

Kit content and related products

| Pneumocystis ELiTe MGB Kit | Pneumocystis ELiTe Standard | Pneumocystis - ELiTe Positive Control |
|---|---|---|
|  X 8 |  X 2 |  X 3 |
| Ready-to-use PCR Mix 8 tubes of 280 µL 96 reactions per kit 7 freeze-thaw cycles | Ready-to-use 4 levels: 10 ⁵ , 10 ⁴ , 10 ³ , 10 ² 2 set of 4 tubes of 200 µL 4 freeze-thaw cycles | Ready-to-use PC 3 tubes of 160 µL 12 reactions per kit 4 freeze-thaw cycles |

Maximum shelf-life: **24 months**

Storage Temperature: **-20 °C**

Other products required not provided in the kit

| | |
|--|--|
| <ul style="list-style-type: none"> ELiTe InGenius instrument: INT030. ELiTe BeGenius instrument: INT040. ELiTe InGenius SP 200: INT032SP200. ELiTe InGenius SP200 Consumable Set: INT032CS. ELiTe InGenius PCR Cassette: INT035PCR. ELiTe InGenius Waste Box: F2102-000. | <ul style="list-style-type: none"> CPE - Internal Control: CTRCPE 300 µL Filter Tips Axigen: TF-350-L-R-S. 1000 µL Filter Tips Tecan: 30180118. |
|--|--|

ELiTe InGenius and ELiTe BeGenius protocol

Table 29

| | | | |
|---------------------------|--------|-----------------------------|-----------|
| › Sample volume | 200 µL | › R/MG PCR Mix volume | 20 µL |
| › CPE volume | 10 µL | › Frequency of controls | 15 days |
| › Total elution volume | 100 µL | › Frequency of calibration | 60 days |
| › Eluate PCR input volume | 20 µL | Unit of quantitative result | Copies/mL |

ELiTe InGenius and ELiTe BeGenius Performances

| Matrix | Limit of Detection | Linearity Range | Diagnostic Sensitivity | Diagnostic Specificity |
|--------|--------------------|-----------------|------------------------|------------------------|
| BAL | 97 cp / mL | 97 – 31,622,777 | 100% 43/43* | 96.6% 56/58* |
| Sputum | 97 cp / mL | 97 – 31,622,777 | 84.2% 32/38* | 97.1% 34/35* |

*confirmed samples/ tested samples

Sample preparation

This product is intended for use on the **ELiTe InGenius** and **ELiTe BeGenius** with the following clinical specimens identified according to laboratory guidelines, and collected, transported, and stored under the following conditions.

| Specimen | Collection requirements | Transport/Storage conditions | | | |
|--|---|---------------------------------|------------|-------------|-------------|
| | | +16 / +26 °C (room temperature) | +2 / +8 °C | -20 ± 10 °C | -70 ± 15 °C |
| Bronchoalveolar Lavage / bronchial aspirate (BAL/BA) | in sterile physiological solution or sterile PBS* | | ≤1 week | ≤ 30 days | ≤ 1 year |
| Sputum (SP) | | | ≤1 week | ≤ 30 days | > 30 days |

*PBS: Phosphate Buffered Saline solution

ELiTe InGenius Procedures

The user is guided step-by-step by the Graphic User Interface (GUI) of ELiTe InGenius software to setup the run. All the steps: extraction, Real-Time PCR and result interpretation are automatically performed. Two operational modes are available: complete run (Extract + PCR) or PCR Only.

Before analysis

| | | |
|--|---|---|
| 1. Switch on ELiTe InGenius. Log in with username and password. Select the mode “ CLOSED ”. | 2. Verify controls: Positive Control and Negative Control in the “Controls” menu. Note: Both must have been run, approved and not expired. | 3. Thaw the PCR Mix and the CTRCPE tubes. Vortex gently. Spin down 5 sec. |
|--|---|---|

Procedure 1 - Complete run: Extract + PCR (e.g., samples)

| | | |
|---|---|--|
| 1. Select “Perform Run” on the touch screen | 2. Verify the extraction volumes: Input: “200 µL”, elution: “100 µL” | 3. Scan the sample barcodes with hand-barcode reader or type the sample ID |
| 4. Select the “Assay Protocol” of interest: PJ ELiTe_BAL_200_100 or PJ ELiTe_SP_200_100 or | 5. Select the method “Extract + PCR” and the sample position: Primary tube or Extraction Tube | 6. Load the PCR Mix and the Internal Control in the Inventory Block |
| 7. Load: PCR Cassette, Extraction cartridge, Elution tube, Tip Cassette, Extraction Tube racks and primary sample racks | 8. Close the door. Start the run | 9. View, approve and store the results |

NOTE

If an Extract Only mode is needed, refer to the instrument user’s manual for procedure.

Procedure 2: PCR Only (e.g., eluates, standards, controls)

| | | |
|--|---|--|
| 1. Select “Perform Run” on the touch screen | 2. Verify the extraction volumes: Input: “200 µL”, elution: “100 µL” | 3. Scan the sample barcodes with hand-barcode reader or type the sample ID |
| 4. Select the “Assay protocol” of interest: PJ ELiTe_PC and PJ ELiTe_NC) | 5. Select the method “PCR Only” and the sample position “Elution Tube” | 6. Load the PCR Mix in the Inventory Block |
| 7. Load: PCR Cassette rack and the Elution tube rack with the extracted nucleic acid | 8. Close the door. Start the run | 9. View, approve and store the results |

ELiTe BeGenius Procedures

The user is guided step-by-step by the Graphic User Interface (GUI) of ELiTe BeGenius software to setup the run. All the steps: extraction, Real-Time PCR and result interpretation are automatically performed. Two operational modes are available: complete run (Extract + PCR) or PCR Only.

Before analysis

| | | |
|--|---|---|
| 1. Switch on ELiTe InGenius. Log in with username and password. Select the mode “ CLOSED ”. | 2. Verify controls: Positive Control and Negative Control in the “Controls” menu. Note: Both must have been run, approved and not expired. | 3. Thaw the PCR Mix and the CTRCPE tubes. Vortex gently. Spin down 5 sec. |
|--|---|---|

Procedure 1 - Complete run: Extract + PCR (e.g., samples)

| | | |
|--|---|---|
| 1. Select "Perform Run" on the touch screen and then click on the run mode «Extract + PCR» | 2. Insert the Sample Rack with the barcoded samples in the Cooler Unit. The barcode scan is already active | 3. Verify the extraction volumes: Input: "200 µL", Eluate: "100 µL" |
| 4. Select the "Assay protocol" of interest (PJ ELITE_Be_BAL_200_100 or PJ ELITE_Be_SP_200_100 or Note: If a second extraction is performed repeat steps from 2 to 4 | 5. Print the labels to barcode the empty elution tubes. Load the tubes in the Elution Rack and insert it in the Cooler Unit | 6. Load the PCR Mix and the Internal Control in the Reagent/Elution Rack and insert it in the Cooler Unit |
| 7. Load "PCR Rack" with "PCR Cassette" and the "Extraction Basket" with the "ELITE InGenius SP 200" extraction cartridges and the required extraction consumables | 8. Close the door. Start the run | 9. View, approve and store the results |

NOTE

If an Extract Only mode is needed, refer to the instrument user's manual for procedure.

Procedure 2: PCR Only (e.g., eluates, controls)

| | | |
|--|--|---|
| 1. Select "Perform Run" on the touch screen and then click on the run mode «PCR Only» | 2. Load the extracted nucleic acid or controls barcoded tubes in the Elution Rack and insert it in the Cooler Unit | 3. Verify the extraction volumes: Input: "200 µL", Eluate: "100 µL" |
| 4. Select the "Assay protocol" of interest (PJ ELITE_Be_PC and PJ ELITE_Be_NC or PJ ELITE_Be_STD) | 5. Load the PCR-Mix in the Reagent/Elution Rack and insert it in the Cooler Unit | 6. Load "PCR Rack" with "PCR Cassette" |
| 7. Close the door. Start the run | 8. View, approve and store the results | |

ELITechGroup S.p.A.
C.so Svizzera, 185, 10149 Torino ITALY
Tel. +39-011 976 191
Fax +39-011 936 76 11
E. mail: emd.support@elitechgroup.com
WEB site: www.elitechgroup.com

