

NOTICE of CHANGE dated 21/12/21

IMPORTANT COMMUNICATION FOR THE USERS OF PRODUCT:

«HSV1 - ELITe Positive Control»







Ref. CTR031PLD

This new revision of the Instruction for Use (IFU) contains the following changes:

- *Update for the use of the product in association with «ELITe BeGenius®» instrument (REF INT040).*

Composition, use and performance of the product remain unchanged.

PLEASE NOTE

| | |
|---|--|
|  | LA REVISIONE DI QUESTO IFU E' COMPATIBILE ANCHE CON LA VERSIONE PRECEDENTE DEL KIT |
|  | THE REVIEW OF THIS IFU IS ALSO COMPATIBLE WITH THE PREVIOUS VERSION OF THE KIT |
|  | CET IFU MIS A JOUR ANNULE ET REMPLACE ET EST PARFAITEMENT COMPATIBLE AVEC LA VERSION PRECEDENTE DU KIT |
|  | LA REVISIÓN DE ESTE IFU ES COMPATIBLE TAMBIÉN CON LA VERSIÓN ANTERIOR DEL KIT |
|  | A REVISÃO DO ESTE IFU ÉTAMBÉM COMPATÍVEL COM A VERSÃO ANTERIOR DO KIT |
|  | DIE REVIEW VON DIESER IFU IST KOMPATIBLE MIT DER VORIGE VERSION VON DEM TEST-KIT |



HSV1 - ELITe Positive Control

plasmid DNA control for qualitative assay

REF CTR031PLD

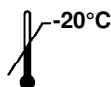


TABLE OF CONTENTS

| | |
|--|--------|
| INTENDED USE | page 1 |
| PRODUCT DESCRIPTION | page 1 |
| MATERIALS PROVIDED IN THE PRODUCT | page 2 |
| MATERIALS REQUIRED BUT NOT PROVIDED IN THE PRODUCT | page 2 |
| OTHER PRODUCTS REQUIRED | page 2 |
| WARNINGS AND PRECAUTIONS | page 2 |
| PROCEDURE | page 3 |
| REFERENCES | page 4 |
| SYMBOLS | page 4 |

INTENDED USE

The «**HSV1 - ELITe Positive Control**» product is intended for use as a positive control in nucleic acids amplification assays for the **detection of the DNA of type 1 Herpes Simplex human virus (HSV1)** with «**HSV1 ELITe MGB® Kit**» product manufactured by ELITechGroup S.p.A.

PRODUCT DESCRIPTION

The product supplies the **Positive Control**, a stabilized solution of plasmid, aliquoted into **two ready-to-use test tubes**. Each test tube contains **160 µL** of solution, sufficient for **4 sessions** in association with the systems «**ELITe InGenius®**» and «**ELITe BeGenius®**», **6 session** in association with the system «**ELITe GALAXY**» and **7 sessions** in association with the other systems validated, considering a manual PCR set-up, mentioned in the instruction for use of the «**HSV1 - ELITe MGB® Kit**» product.

Plasmid contains a region of the gene codifying the **glycoprotein D (gpD)** of HSV1. This plasmid is detected by the real time amplification reaction, attesting the reaction ability to detect the HSV1 DNA.

The product is sufficient for **8 separate analytic sessions** in association with the systems «**ELITe InGenius®**» and «**ELITe BeGenius®**», **12 separate analytic sessions** in association with the system «**ELITe GALAXY**» and **14 separate analytic sessions** in association with the other systems considering a manual PCR set-up, by using 20 µL for reaction.

HSV1 - ELITe Positive Control
plasmid DNA control for qualitative assay

REF CTR031PLD

MATERIALS PROVIDED IN THE PRODUCT

| Component | Description | Quantity | Classification and Labelling |
|-------------------------|------------------|------------|------------------------------|
| HSV1 - Positive Control | Plasmid solution | 2 x 160 µL | - |

MATERIALS REQUIRED BUT NOT PROVIDED IN THE PRODUCT

- Laminar airflow hood.
- Disposable powderless gloves in nitrile or similar material.
- Vortex mixer.
- Bench microcentrifuge (12,000 - 14,000 RPM).
- Micropipettes and sterile tips with aerosol filter or sterile positive displacement tips (2-20 µL, 5-50 µL, 50-200 µL).
- Molecular biology grade water.
- Programmable thermostat with optical fluorescence detection system 7300 Real Time PCR System or 7500 Fast Dx Real-Time PCR Instrument calibrated following manufacturer's instructions.
- Programmable thermostat with optical fluorescence detection system cobas z 480 analyzer, calibrated following manufacturer's instructions.

OTHER PRODUCTS REQUIRED

The reagents for real time amplification and the consumables **are not** included in this product.

To perform real time amplification, the product «**HSV1 ELITe MGB® Kit**» (ELITechGroup S.p.A, ref. RTS031PLD), it is required. This is a complete and ready for use reaction mixture for real time amplification in a stabilising solution.

In association with «**ELITe InGenius®**» (ELITechGroup S.p.A, ref. INT030) and «**ELITe BeGenius®**» (ELITechGroup S.p.A, ref. INT040) it is required the use of generic product «**ELITe InGenius® PCR Cassette**» (ELITechGroup S.p.A, ref. INT035PCR). These are dedicated consumables for Real Time PCR reactions.

When a 7300 Real-Time PCR System is used, it is required the use of generic product «**Q - PCR Microplates**» (ELITechGroup S.p.A., ref. RTSACC01), microplates with 0.2 mL wells and adhesive sealing sheets for real time amplification.

When a 7500 Fast Dx Real-Time PCR Instrument is used, it is required the use of generic product: «**Q - PCR Microplates Fast**» (ELITechGroup S.p.A., ref. RTSACC02), microplates with 0.1 mL wells and adhesive sealing sheets for real time amplification.

When a cobas z 480 analyzer is used, it is required the use of generic product «**AD-plate 0.3ml**» (Roche, ref. 05232724001), microplates with 0.3 mL wells and adhesive sealing sheets for real time amplification.

WARNINGS AND PRECAUTIONS

This product is exclusively for *in vitro* use.

Warnings and general precautions

Handle and dispose of all biological samples as if they were able to transmit infective agents. Avoid direct contact with the biological samples. Avoid splashing or spraying. The materials that come into contact with the biological samples must be treated for at least 30 minutes with 3% sodium hypochlorite or autoclaved for one hour at 121°C before disposal.

Handle and dispose of all reagents and all materials used to carry out the assay as if they were able to transmit infective agents. Avoid direct contact with the reagents. Avoid splashing or spraying. Waste must be handled and disposed of in compliance with adequate safety standards. Disposable combustible material must be incinerated. Liquid waste containing acids or bases must be neutralised before disposal.

HSV1 - ELiTe Positive Control
plasmid DNA control for qualitative assay

REF CTR031PLD

Wear suitable protective clothes and gloves and protect eyes and face.
Never pipette solutions by mouth.
Do not eat, drink, smoke or apply cosmetic products in the work areas.
Carefully wash hands after handling samples and reagents.
Dispose of leftover reagents and waste in compliance with the regulations in force.
Carefully read all the instructions provided in the product before running the assay.
While running the assay, follow the instructions provided in the product.
Do not use the product after the indicated expiry date.
Only use the reagents provided in the product and those recommended by the manufacturer.
Do not use reagents from different batches.
Do not use reagents from other manufacturers.

Warnings and precautions for molecular biology

Molecular biology procedures, such as nucleic acids extraction, amplification and detection, require qualified and trained staff to avoid the risk of erroneous results, especially due to the degradation of nucleic acids contained in the samples or sample contamination by amplification products.

When amplification session is manually setup, it is necessary to have available separate areas for the extraction / preparation of amplification reactions and for the amplification / detection of amplification products. Never introduce an amplification product in the area designated for extraction / preparation of amplification reactions.

When amplification session is manually setup, it is necessary to have available lab coats, gloves and tools which are exclusively used for the extraction / preparation of the amplification reactions and for the amplification / detection of amplification products. Never transfer lab coats, gloves or tools from the area designated for the amplification / detection of amplification products to the area designated for the extraction / preparation of the amplification reactions.

The samples must be exclusively used for this type of analysis. Samples must be handled under a laminar airflow hood. Tubes containing different samples must never be opened at the same time. Pipettes used to handle samples must be exclusively used for this specific purpose. The pipettes must be of the positive displacement type or be used with aerosol filter tips. The tips used must be sterile, free from DNases and RNases, free from DNA and RNA.

The reagents must be handled under a laminar airflow hood. The reagents required for amplification must be prepared in such a way that they can be used in a single session. The pipettes used to handle the reagents must be exclusively used for this purpose. The pipettes must be of the positive displacement type or be used with aerosol filter tips. The tips used must be sterile, free from DNases and RNases, free from DNA and RNA.

Amplification products must be handled in such a way as to reduce as much as possible dispersion into the environment in order to avoid the possibility of contamination. The pipettes used to handle amplification products must be exclusively used for this purpose.

Warnings and precautions specific for the components

The **Positive Control** can be frozen and thawed for no more than **seven times**. Further freezing / thawing cycles may cause a loss in titre.

The **Positive Control** can be left on board on the «**ELiTe InGenius®**» and on the «**ELiTe BeGenius®**» instruments for up to **four work sessions of three hours each** ("Extract + PCR" mode).

PROCEDURE

The «**HSV1 - ELiTe Positive Control**» product must be used with the complete reaction mixture of the «**HSV1 ELiTe MGB® Kit**» product.

Before use, take and thaw the **HSV1 - Positive Control** tubes. Mix gently, spin down the contents for 5 seconds and keep them on ice.

The **HSV1 - Positive Control** is ready to use: **20 µL** must be directly added to the reaction mixture.

The complete procedure, the performance characteristics and procedure limitations of the complete assay are described in detail in the instructions for use of the «**HSV1 ELiTe MGB® Kit**» product.

HSV1 - ELiTe Positive Control
plasmid DNA control for qualitative assay

REF CTR031PLD






Note: The results of the Positive Control amplification will be stored by the instrument «**ELiTe InGenius®**» and «**ELiTe BeGenius®**» and used to create a control chart. For each batch of product «**HSV1 ELiTe MGB® Kit**» the amplification of Positive Control is required. The stored results of the Positive Control amplification will expire after **15 days**.

Note: The **Positive Control** can be frozen and thawed for no more than **seven times**. The **Positive Control** can be left on board on the «**ELiTe InGenius®**» and «**ELiTe BeGenius®**» instruments for up to **four work sessions of three hours each** ("Extract + PCR" mode).

REFERENCES

E. Aurelius et al. (1993) *J. Med. Virology* **39**: 179 - 186

SYMBOLS

| | |
|---|---|
| REF | Catalogue Number. |
|  | Upper limit of temperature. |
| LOT | Batch code. |
|  | Use by (last day of month). |
| IVD | <i>in vitro</i> diagnostic medical device. |
| CE | Fulfilling the requirements of the European Directive 98/79/EC for <i>in vitro</i> diagnostic medical device. |
|  | Contains sufficient for "N" tests. |
|  | Attention, consult instructions for use. |
| CONT | Contents. |
|  | Manufacturer. |

ELiTe MGB®, the ELiTe MGB® logo and ELiTe InGenius® and ELiTe BeGenius® are registered as trademarks within the European Union.