

# RF CALIBRATOR

REF IRFA-0042



PITV-IRFACal-v4 (06/2021)

146 IU/mL

LOT


22-4015



2024-08


Traceability : NIBSC 64/002

## Preparation

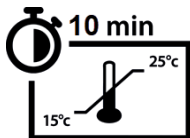


RF Cal	μL	-	25	50	100	200	400
NaCl 9 g/L	μL	400	375	350	300	200	-
Factor		0	1/16	1/8	1/4	1/2	1

## Precautions

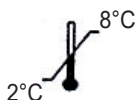

[www.elitechgroup.com](http://www.elitechgroup.com)

selectra  
**PRO**

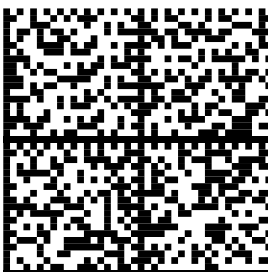

RF CALIBRATOR 0  
S80 PITV-IRFACal  
22-4015 2024-08



ELITech Clinical Systems SAS - Zone Industrielle - 61500 SEES FRANCE